

**THE HASHEMITE UNIVERSITY**

**DEPARTMENT OF CHEMISTRY**



LABORATORY MANUAL  
FOR  
EXPERIMENTAL CHEMISTRY FOR  
ENGINEERING  
CHEM. 110103108

**Version 6**

**(2020)**

**Prepared by:**

***Dr. Khalid Al-Ani***  
***Dr. Adnan Abu-Surrah***

***Dr. Hamzeh Abdel-Halim***  
***Dr. Nayif Masoud***



**Chemistry Laboratory for Engineering Students  
(110103108)**

**Course Content :**

<b><u>Exp. No</u></b>	<b><u>Experiment Title</u></b>	<b><i>Page</i></b>
1.	Safety and Equipment	3
2.	Empirical Formula of a Compound	11
3.	Limiting Reactant	18
4.	Identification of a Compound: Physical Properties	25
5.	Tests of Cations and Anions	32
6.	Molar Mass of a Volatile Liquid	37
7.	Determination of the Molar Volume of Hydrogen Gas	43
8.	Colligative Properties: Molar Mass Determination	48
9.	Calorimetry	56
10.	Electrochemistry	61
11.	Chemical Thermodynamics	66
12.	pH-Hydrolysis	70
13.	Acid-Base Titration	75

1 1A <b>H</b> 1.008	2 2A <b>He</b> 4.003	3 <b>Li</b> 6.941	4 <b>Be</b> 9.012	5	6	7	8	9	10	11	12	13 3A <b>B</b> 10.81	14 4A <b>C</b> 12.01	15 5A <b>N</b> 14.01	16 6A <b>O</b> 16.00	17 7A <b>F</b> 19.00	18 8A <b>Ne</b> 20.18
11 <b>Na</b> 22.99	12 <b>Mg</b> 24.31											13 <b>Al</b> 26.98	14 <b>Si</b> 28.09	15 <b>P</b> 30.97	16 <b>S</b> 32.07	17 <b>Cl</b> 35.45	18 <b>Ar</b> 39.95
19 <b>K</b> 39.10	20 <b>Ca</b> 40.08	21 <b>Sc</b> 44.96	22 <b>Ti</b> 47.88	23 <b>V</b> 50.94	24 <b>Cr</b> 52.00	25 <b>Mn</b> 54.94	26 <b>Fe</b> 55.85	27 <b>Co</b> 58.93	28 <b>Ni</b> 58.69	29 <b>Cu</b> 63.55	30 <b>Zn</b> 65.38	31 <b>Ga</b> 69.72	32 <b>Ge</b> 72.59	33 <b>As</b> 74.92	34 <b>Se</b> 78.96	35 <b>Br</b> 79.90	36 <b>Kr</b> 83.80
37 <b>Rb</b> 85.47	38 <b>Sr</b> 87.62	39 <b>Y</b> 88.91	40 <b>Zr</b> 91.22	41 <b>Nb</b> 92.91	42 <b>Mo</b> 95.94	43 <b>Tc</b> (98)	44 <b>Ru</b> 101.1	45 <b>Rh</b> 102.9	46 <b>Pd</b> 106.4	47 <b>Ag</b> 107.9	48 <b>Cd</b> 112.4	49 <b>In</b> 114.8	50 <b>Sn</b> 118.7	51 <b>Sb</b> 121.8	52 <b>Te</b> 127.6	53 <b>I</b> 126.9	54 <b>Xe</b> 131.3
55 <b>Cs</b> 132.9	56 <b>Ba</b> 137.3	57 <b>La*</b> 138.9	72 <b>Hf</b> 178.5	73 <b>Ta</b> 180.9	74 <b>W</b> 183.9	75 <b>Re</b> 186.2	76 <b>Os</b> 190.2	77 <b>Ir</b> 192.2	78 <b>Pt</b> 195.1	79 <b>Au</b> 197.0	80 <b>Hg</b> 200.6	81 <b>Tl</b> 204.4	82 <b>Pb</b> 207.2	83 <b>Bi</b> 209.0	84 <b>Po</b> (209)	85 <b>At</b> (210)	86 <b>Rn</b> (222)
87 <b>Fr</b> (223)	88 <b>Ra</b> 226	89 <b>Ac<sup>†</sup></b> (227)															

*Lanthanides	58 <b>Ce</b> 140.1	59 <b>Pr</b> 140.9	60 <b>Nd</b> 144.2	61 <b>Pm</b> (145)	62 <b>Sm</b> 150.4	63 <b>Eu</b> 152.0	64 <b>Gd</b> 157.3	65 <b>Tb</b> 158.9	66 <b>Dy</b> 162.5	67 <b>Ho</b> 164.9	68 <b>Er</b> 167.3	69 <b>Tm</b> 168.9	70 <b>Yb</b> 173.0	71 <b>Lu</b> 175.0
†Actinides	90 <b>Th</b> 232.0	91 <b>Pa</b> (231)	92 <b>U</b> 238.0	93 <b>Np</b> (237)	94 <b>Pu</b> (244)	95 <b>Am</b> (243)	96 <b>Cm</b> (247)	97 <b>Bk</b> (247)	98 <b>Cf</b> (251)	99 <b>Es</b> (252)	100 <b>Fm</b> (257)	101 <b>Md</b> (258)	102 <b>No</b> (259)	103 <b>Lr</b> (260)

## **Experiment (1)**

### **Safety and Equipments**

Welcome to the Chemistry lab! Chemistry is an experimental science and you will learn many important things by doing experiments. In this first laboratory, we want you to get familiar with safety rules in the laboratory. Also, you will be introduced to the laboratory so that you'll get familiar with the lab and your equipment. **There is no grade for this first lab session!**

#### **OBJECTIVES OF THE PRACTICAL COURSE IN CHEMISTRY**

1. To learn basic laboratory skills with standard chemical equipment.
2. To carry out experiments involving laboratory chemicals in a safe and correct way with minimum waste.
3. To learn to observe and report the results of simple chemical experiments correctly.
4. To learn to explain observations using the concepts you have learned in the lecture course.

**Read the following safety rules before going into the laboratory and abide by them.**

**Wear safety glasses at all times. Prescription glasses are adequate but contact lenses are not adequate.**

- 1- Always wear your laboratory coat. Do not wear clothing that hinders free movement of your hands or hangs loose outside your laboratory coat.
- 2- Do not work in a laboratory if no lecturer or technician is present. Read the experimental instructions carefully before starting the work. Especially note any precautions that must be taken.
- 3- Never eat, drink, or smoke in the laboratory. Never taste chemicals. Wash your hands well before leaving the laboratory. Also, wash your hands or any part of your body immediately with water when it comes in contact with chemicals.
- 4- Do not use your mouth to fill a pipette. There are special bulbs for this purpose.
- 5- Note the position of safety equipment like fire extinguishers, eye washers, and first aid boxes. Report all accidents immediately to a staff member or technician.
- 6- Use the fume hood when handling strong-smelling or irritating chemicals.
- 7- Be careful about discarding away wastes. Always follow instructions.

Do not dispose of solids into the sink. Do not leave glassware or any other solid materials, including filter papers, in the sink. Put broken glassware into the labeled buckets.


Some waste liquids must be stored into special bottles, not disposed of in the sink. A staff member or technician will help you.

- 8- Do not leave a lit burner unattended. Always stay clear form the flame.
- 9- When heating anything in a test tube, do not point the mouth of the test tube towards yourself or towards any other person.
- 10- Before leaving the laboratory turn off any water taps and burners and dispose of solid waste in the correct container. Also, wash all apparatus used and clean up the bench top.
- 11- Keep your bench clean and tidy while you are working. Clean up any spills or broken glass immediately. Keep your books and papers away from water, chemicals, and flames. Position your apparatus on the bench so that it is convenient and comfortable to use. Keep unused equipment out of the way, so that you do not knock it over.
- 12- If you are in any doubt about anything, ask the staff member or technician for advice.
- 13- Request the assistance of your instructor/technician if and when you suffer a cut or a burn **or face a dangerous** situation.


## SOME LABORATORY APPARATUS




Beaker

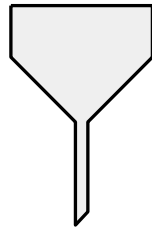


Büchner funnel

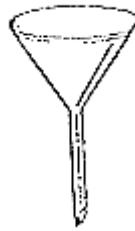


Funnel

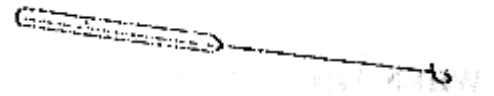

Glass rod with platinum wire  
(for flame testing)

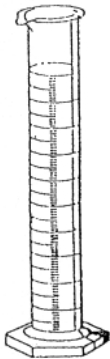


Graduated cylinder

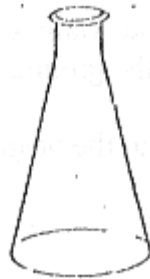


Erlenmeyer flask



Dropping bottle

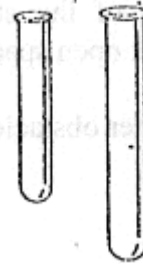

Test tubes

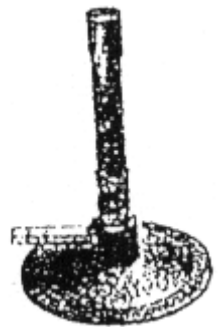

Bunsen burner



Stirring rod



Dropper



Safety glasses  
(Goggles)




Test tube rack



Plastic wash bottle



Clamp holder

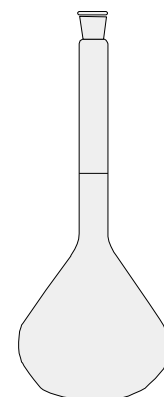


Volumetric flask



Burette




Tripod

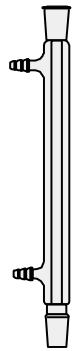

Condenser



Stand

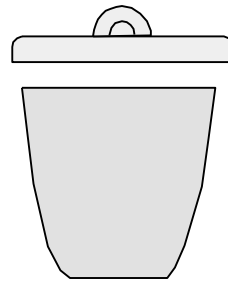


Crucible




Volumetric pipette




Reagent bottle

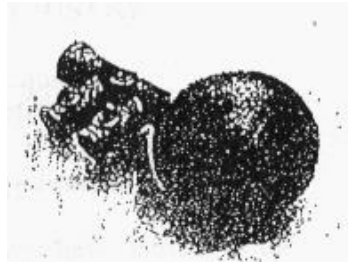

Pipette filler

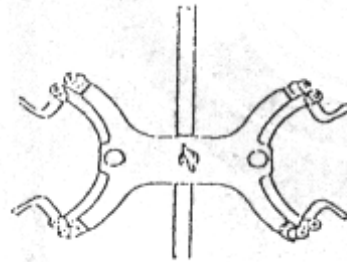

Burette clamp

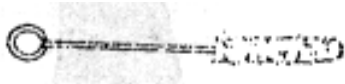

Test tube brush



Metal spatula



Test tube holder




Ring clamp



Ceramic square



Forceps

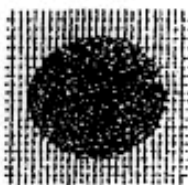

Wire gauze



Watch glass



Crucible tongs



**Top-loading  
balance**



## FIRE REGULATIONS

Fire is a serious hazard in any laboratory. Always assume that an alarm is real. There is no need to leave the building, but only when you are instructed by your instructor to do so.

### WHEN THE CONTINUOUS ALARMS SOUND:

1. Do not stay to collect personal belongings.
2. Follow the instructions of the staff and walk outside through the nearest exit to the nearest open space on the ground level.
3. If you meet smoke or other obstacles, go to the nearest clear exit.

### WHEN THE INTERMITTENT ALARMS SOUND:

There is no need to leave the building, but be ready to leave if the alarm becomes continuous.

## INSTRUCTIONS FOR ALL EXPERIMENTS

You may lose marks if you do not follow these instructions:


1. Wear a laboratory coat and safety glasses in the laboratory at all times.
2. Write all your results in the spaces provided in the laboratory manual at the time of the experiment.
3. Write all your results in blue or black ink, not in pencil.
4. If you make a mistake and you need to change an entry, don't erase it, overwrite it, or use correction fluid. Cross out the old value and write the new value as follows:

$$\text{Mass of water} = \begin{array}{r} 10.15 \text{ g} \\ \times \\ 10.51 \text{ g} \end{array}$$

5. Complete your experimental report, including calculations, in the laboratory. You should bring a calculator to the laboratory.
6. **Submit your report to your instructor for marking before you clean up your bench** and leave the laboratory.
7. When you have finished the experiment, clean your apparatus, and return everything that you took from your locker. **Leave your bench clean and tidy** when you go.
8. You may not work in the laboratory after the end of the laboratory period. Leave the laboratory on time so as to be able to reach your next class on time.

On this drawing, please mark the locations of the following equipment. There may be more than one item of each kind. Please mark them all. Walk around the lab to find each item. **You may need help from your instructor** or a technician. Here is the list of items to find.

- |                       |                                    |
|-----------------------|------------------------------------|
| 1. All exit doors     | 5. Fire blankets                   |
| 2. Fume hoods         | 6. Emergency shower                |
| 3. Fire extinguishers | 7. First aid kits                  |
| 4. Eye wash stations  | 8. Poster with safety instructions |



### A. Safety

A chemistry lab could be a dangerous place to work in. To make sure that nobody gets hurt we ask you to follow a few simple instructions. You should already have read the safety instructions on pages 2,3 and 7. Please indicate that you have done so by signing your name below:

I have read and understood the 'Safety in the Laboratory' and the 'Fire Regulations' on pages 2,3 and 7. I agree to follow these rules.

Name \_\_\_\_\_ ID No. \_\_\_\_\_

The lab has many safety features that can be used in an emergency. On the previous page you will find a drawing of the lab. Please mark the locations of all the safety equipment listed there.

## **B. Equipment**

As in all new fields of study, there are some new words and equipment to get to know in the chemistry lab. Please take a few minutes to get familiar with the glassware and equipment you will be using during this course. Here is what to do:

Start by taking all items from your drawer and your locker and place them on the bench. Now take each item one by one and see if you can find its name on the list on pages 4, 5 and 6. If you find the item on the list, circle the drawing. Just ask if you need help.

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## **Experiment (2)**

### **Empirical Formula of a Compound**

#### **PreLab**

1. How many grams of magnesium combine with 1.50 g of chloride ions in  $\text{MgCl}_2$ ?

---

---

---

---

---

2. If 11.80 g of iron reacts with 5.06 g of oxygen. Determine the empirical formula of the resulting oxide?

---

---

---

---

---

---

---

---

---

---

## Experiment (2)

### Empirical Formula of a Compound

#### Objective:

- To determine the Chemical formula of magnesium oxide.

#### Introduction:

The empirical formula of a compound is the smallest whole – number ratio of atoms in the compound. The following are examples of empirical formulas: Fe<sub>2</sub>O<sub>3</sub>, ZnO, Mg<sub>3</sub>N<sub>2</sub>, CH<sub>3</sub>, CH<sub>2</sub>, CH, CO<sub>2</sub>, MgS, and MnO<sub>2</sub>. The empirical formula of a compound is usually determined either by chemical analysis (decompositions) or by synthesis. In the chemical analysis method, a known mass of the compound is decomposed (or burned) and the masses of the products are determined. These masses are converted to number of moles of the elements present in the compound. From the number of moles, the relative number of atoms of each element is determined, and this will give the empirical formula.

#### Example 1:

Nicotine is a compound containing C, H, and N. A 2.50 g sample of the compound is burned and produces 6.78 g of CO<sub>2</sub>, 1.94 g of H<sub>2</sub>O, and 0.43 g of N<sub>2</sub>. what is the empirical formula of nicotine?

First we calculate the number of moles of each element present.

$$\text{mole C} = \text{mole CO}_2 = \frac{6.78\text{g}}{44.0\text{g/mol}} = 0.154 \text{ mole}$$

$$\text{mole H} = 2 \times \text{mole H}_2\text{O} = 2 \times \frac{1.94\text{g}}{18.0\text{g/mol}} = 0.216 \text{ mole}$$

$$\text{mole N} = 2 \times \text{mole N}_2 = 2 \times \frac{0.43\text{g}}{28.0\text{g/mol}} = 0.031 \text{ mole}$$

ratio of moles = ratio of atoms

The empirical formula is C<sub>0.154</sub> : H<sub>0.216</sub> : N<sub>0.031</sub>

Dividing by 0.031, (the smallest number) this will give the smallest whole-number ratio:

$C_5H_7N$  = the empirical formula

In the chemical synthesis method, a known mass of a metal is combined with a non metal. The product is weighed and the mass of nonmetal is determined. From the mass of the metal and nonmetal, the number of moles is calculated and then the empirical formula is determined.

**Example 2:**

When 0.288 g of P is burned, 0.660 g of white phosphorus oxide is obtained, determine the empirical formula of this oxide?

$$\text{mass of oxygen} = 0.660 - 0.188 = 0.372 \text{ g}$$

$$\text{mole P} = \frac{0.288 \text{ g}}{31.0 \text{ g/mol}} = 0.0093 \text{ mole}$$

$$\text{mole O} = \frac{0.372 \text{ g}}{16.0 \text{ g/mol}} = 0.0233 \text{ mole}$$

$$\text{mole P: mole O} = 0.0093: 0.0233 = 2 : 5$$

Empirical formula is  $P_2O_5$

In today's experiment you will prepare magnesium oxide by burning a known mass of magnesium in air ( $O_2$ ). From the mass of product, and magnesium, the empirical formula of magnesium oxide is determined.

**Apparatus and Chemicals:**

Balance, Bunsen burner, Crucible and cover, Clay triangle, Crucible tongs, Iron ring, clamps, stand, Mg ribbon and Steel wool.

## Experimental Procedure:

Record all observations and data directly on report sheet **IN INK**.

1. Place a clean, dry, partially covered crucible on the clay triangle. Heat with hot flame for two or three minutes. This step removes any moisture from the crucible. See fig (1) for correct assembly.
2. Turn off the flame and allow the crucible and cover to cool on the clay triangle.
3. **WHEN THE CRUCIBLE STOPS FEELING WARM TO YOUR HAND**, bring the crucible and cover to a balance and weigh them to two decimal places ( $\pm 0.01$  g). If the crucible is too hot you will not get an accurate weight.
4. Polish about 2 cm of Mg ribbon and cut it into small pieces. Use steel wool for polishing. Place about 0.25 g of the element in the crucible and weigh the crucible + Mg accurately to ( $\pm 0.01$  g).
5. Replace cover and heat the crucible and contents slowly at the beginning. Lift the cover occasionally to introduce oxygen to react with Mg. Do not lift the cover widely, because too much air will burn the Mg brightly (with flame). If this happens cover the crucible immediately.
6. Continue heating until all Mg is converted to ash, then remove the cover and heat the crucible until the crucible bottom becomes red.
7. Remove the burner and Cool the crucible, add a few drops of distilled water to decompose any magnesium nitride,  $\text{Mg}_3\text{N}_2$ , that may have formed during combustion.
8. Heat the open crucible slowly until the ash is completely dry.
9. Allow the crucible to cool for 2 minutes on the clay triangle. Then place the covered crucible in the desiccators and cool to **room temperature**, about 15 to 20 minutes. **DO NOT PLACE THE HOT CRUCIBLE ON THE LAB BENCH AS IT WILL BURN THE BENCH.**
10. When the crucible cools to room temperature (that is, not warm to the touch), transfer it with the cover to the balance and weigh them to two decimal places.

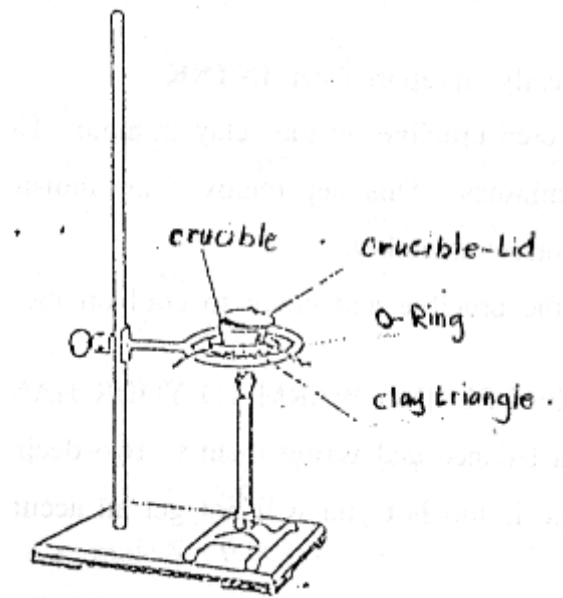


Fig. (1): Setup for magnesium oxide determination.



Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (2)

### Empirical Formula of a Compound

#### Report Sheet

##### Data:

	Trial I	Trial II
1. Mass of empty crucible	_____	_____
2. Mass of Mg	_____	_____
3. Mass of Crucible and oxide	_____	_____
4. Mass of magnesium oxide	_____	_____
5. Mass of oxygen	_____	_____
6. Moles of Mg (FM of Mg = 24.3 g /mole)	_____	_____
7. Moles of oxygen (FM of O = 16.0 g/ mole)	_____	_____
8. Mole Mg: Mole O	_____	_____
9. Empirical formula of magnesium oxide	_____	_____
10. % of O in the compound	_____	_____
11. % of Mg in the compound	_____	_____

**Show your calculation for steps of 6-11 :-**

**Experiment (2)**  
**Empirical Formula of a Compound**

**PostLab**

1. What is the effect of  $Mg_3N_2$  if it is not decomposed on the reported (**Mg to oxygen**) mole ratio?

---

2. Describe the effect of each of the following factors, whether increases, decreases, or has no effect on the reported value of (Mg to O) mole ratios:

a. If carbon is deposited on the crucible's surface (because of improper heating) and the crucible with contents is weighed without removing the carbon residue by further heating.

---

b. If the magnesium oxide ash is not dried completely.

---

c. If rapid oxidation of magnesium occurs by too much air, and some Mg is lost.

---

d. If air is not sufficient to react with all the magnesium.

---

3. Explain how magnesium nitride is formed during heating? Write a chemical equation for its formation.

---

4. Explain, by a chemical equation, the effect of the added water in the decomposition of magnesium nitride.

---

---


Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

### **Experiment (3)**

### **Limiting Reactant**

### **Prelab**

1. A 25.00 g sample of  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$  reacts with excess  $\text{BaCl}_2 \cdot \text{H}_2\text{O}$ . If the mass of  $\text{Ba}_3(\text{PO}_4)_2$  obtained is 17.56 g. Calculate the % yield of  $\text{Ba}_3(\text{PO}_4)_2$ ?

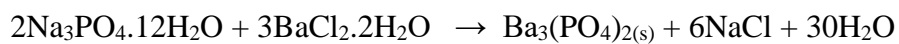

---

---

---

---

---

---

---

---

## Experiment (3)

### Limiting Reactant

#### Objectives:

- To determine the limiting reactant in a salt mixture.
- To determine the percent composition of a salt mixture.

#### Introduction:

In a chemical reaction, the reacting materials are often mixed in stoichiometric quantities (according to fixed mole ratios), i.e., in exactly the correct amounts so that all reactants are used up during the reaction.

When reactants are not mixed in the correct mole ratio, only limited amount of product is formed, and this amount is determined by the reactant present in less amount and thus is consumed first in a reaction. This reactant is called the limiting reactant.


**Limiting reactant:** It is the reactant that is consumed first, and thus determines the amount of product formed.

**The theoretical yield** of a product is the maximum amount that can be produced with a given amount of the limiting reactant. The **actual yield** is often less than the theoretical yield. The % yield is represented as:

$$\% \text{ yield} = \frac{\text{Actual yield}}{\text{Theoretical yield}} \times 100\%$$

Let us discuss the concept of limiting reactant by examining the reaction between  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$  and  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$ .


This is a precipitation reaction which produces  $\text{Ba}_3(\text{PO}_4)_2(\text{s})$ . The equation in molecular form is:

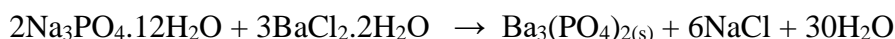

and the ionic equation is:



The net ionic equation of the reaction is:

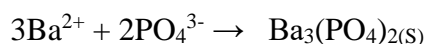

In this reaction:

3 moles of  $\text{Ba}^{2+}$  ions react with 2 moles of  $\text{PO}_4^{3-}$  ions, or


3 moles of  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$  react with 2 moles of  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$

### Example 1:

A mixture containing 40.00 g of  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$  and 30.00 g of  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$  is dissolved in water. A precipitate of  $\text{Ba}_3(\text{PO}_4)_2$  weighing 22.65 g is produced. Calculate the % yield of  $\text{Ba}_3(\text{PO}_4)_2(\text{s})$ .

(FM) of  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O} = 380.2 \text{ g/mole}$ ; (FM) of  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O} = 244.2 \text{ g/mole}$ .

The net ionic equation is:

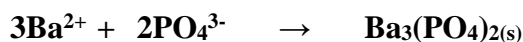

mole  $\text{PO}_4^{3-} = \text{mole } \text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O} = 40.00\text{g} / 380.2 \text{ g mole}^{-1} = 0.105 \text{ mole}$

mole  $\text{Ba}^{2+} = \text{mole } \text{BaCl}_2 \cdot 2\text{H}_2\text{O} = 30.00 \text{ g} / 244.2 \text{ g mole}^{-1} = 0.123 \text{ mole}$

0.105 mole  $\text{PO}_4^{3-}$  requires  $0.105 \times 3/2 = 0.158 \text{ mole } \text{Ba}^{2+}$ . Since we have only 0.123 mole of  $\text{Ba}^{2+}$ , so the **limiting reactant** is  $\text{Ba}^{2+}$ , and the **reactant in excess** is  $\text{PO}_4^{3-}$ .

The **theoretical yield** of product is calculated from the limiting reactant as follows:

$$\text{mole } \text{Ba}_3(\text{PO}_4)_2 = \text{moles } \text{Ba}^{2+} \times \frac{1 \text{ mole } \text{Ba}_3(\text{PO}_4)_2}{3 \text{ mole } \text{Ba}^{2+}}$$

$$= 0.123 \times \frac{1}{3}$$


$$= 0.041 \text{ mole}$$

mass of  $\text{Ba}_3(\text{PO}_4)_2 = 0.041 \text{ mole} \times 602.2 \text{ g mole}^{-1} = 24.69 \text{ g}$  (this is the theoretical yield)

$$\% \text{ yield} = \frac{22.65 \text{ g}}{24.69 \text{ g}} \times 100 = 91.7 \%$$

### Example 2:

10.0 g of a unknown mixture containing  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$  and  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$  is dissolved in distilled water. The mass of  $\text{Ba}_3(\text{PO}_4)_2$  precipitated is 3.50 g. Calculate the % of each salt present in the mixture if  $\text{BaCl}_2$  is the limiting reactant?

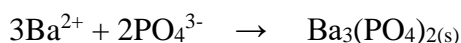

$$\text{mole } \text{Ba}_3(\text{PO}_4)_2 = \frac{3.50 \text{ g}}{602.2 \text{ g/mol}} = 5.81 \times 10^{-3} \text{ mole}$$

$$\text{mole } \text{Ba}^{2+} = \text{mole } \text{BaCl}_2 \cdot 2\text{H}_2\text{O} = 3 \times (\text{mole } \text{Ba}_3(\text{PO}_4)_2)$$

$$\text{mole } \text{BaCl}_2 \cdot 2\text{H}_2\text{O} = 5.81 \times 10^{-3} \times 3 = 1.74 \times 10^{-2} \text{ mole}$$

$$\text{mass of } \text{BaCl}_2 \cdot 2\text{H}_2\text{O} = 1.74 \times 10^{-2} \text{ mole} \times 244.2 \text{ g mole}^{-1} = 4.25 \text{ g}$$

$$\text{mass of } \text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O} = 10.0 - 4.25 = 5.75 \text{ g}$$

$$\% \text{ BaCl}_2 \cdot 2\text{H}_2\text{O} = \frac{4.25}{10.0} \times 100 = 42.5\%$$

$$\% \text{ Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O} = \frac{5.75}{10.0} \times 100 = 57.5\%$$

In the following experiment you will receive an unknown mixture of  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$  and  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$ . The mixture is weighed and then dissolved in water, the precipitate formed is dried and weighed. From the weights of mixture and precipitate, the composition of the unknown mixture is determined.

### Apparatus and Chemicals:

Balance	400- mL beaker
Bunsen Burner	100-mL beakers
Iron ring	Filter flask
Ring stand and clamps	Buchner funnel
Wire gauze	Filter paper (Whatman No. 42)
Watch glass	$\text{BaCl}_2$
Stirring rod	$\text{Na}_3\text{PO}_4$
Rubber policeman	Unknown mixture of ( $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$ ) and ( $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$ )

### Experimental Procedure:

**Two trials** are required. Weigh **duplicate** unknown salt mixtures and perform the following steps for each.

#### (1) Precipitation of $\text{Ba}_3(\text{PO}_4)_2$ :

1. Weigh on a weighing paper about 1.0 g of an **unknown mixture** containing  $\text{BaCl}_2 \cdot 2\text{H}_2\text{O}$  and  $\text{Na}_3\text{PO}_4 \cdot 12\text{H}_2\text{O}$ . Report the mass to  $\pm 0.01$  g.
2. Transfer to a 400 mL beaker, add 200 mL of distilled water and **stir the mixture for one minute**.
3. Cover the beaker with a watch glass and **heat gently** with stirring for 30 minute. **Do not boil the solution**, keep the temperature at about 80-90°C
4. Allow the precipitate to settle, then decant two 50mL volume of the supernatant liquid into two 100 mL beakers, labeled beaker 1 & 2 to be used in part (II).
5. Weigh a filter paper ( $\pm 0.01$  g) and seal it into the Buchner funnel with a small amount of distilled water. While the solution is hot, **filter the precipitate** by vacuum as shown in Fig.1. Transfer any precipitate remaining in the beaker into the funnel with a policeman and hot water.
6. **Wash the precipitate with two portions of distilled water**, remove the filter paper and put it in the oven and dry the precipitate at a temperature of 100°C for about 1 hour.
7. Remove the filter paper from the oven and weight the precipitate and paper to ( $\pm 0.01$  g)

**(2) Determination of the limiting reactant:**

The limiting reactant is determined using the solutions in beakers 1 & 2 as follows:

**Test for excess  $\text{PO}_4^{3-}$ :** add 2 drops of 0.5 M  $\text{BaCl}_2$  to the solution in beaker 1. If a precipitate is formed, then  $\text{PO}_4^{3-}$  is in excess and  $\text{Ba}^{2+}$  is the limiting reactant.

**Test for excess  $\text{Ba}^{2+}$ :** add 2 drops of 0.5 M  $\text{Na}_3\text{PO}_4$  to the solution in beaker 2. If a ppt. is formed, then  $\text{Ba}^{2+}$  is in excess and  $\text{PO}_4^{3-}$  is the limiting reactant.

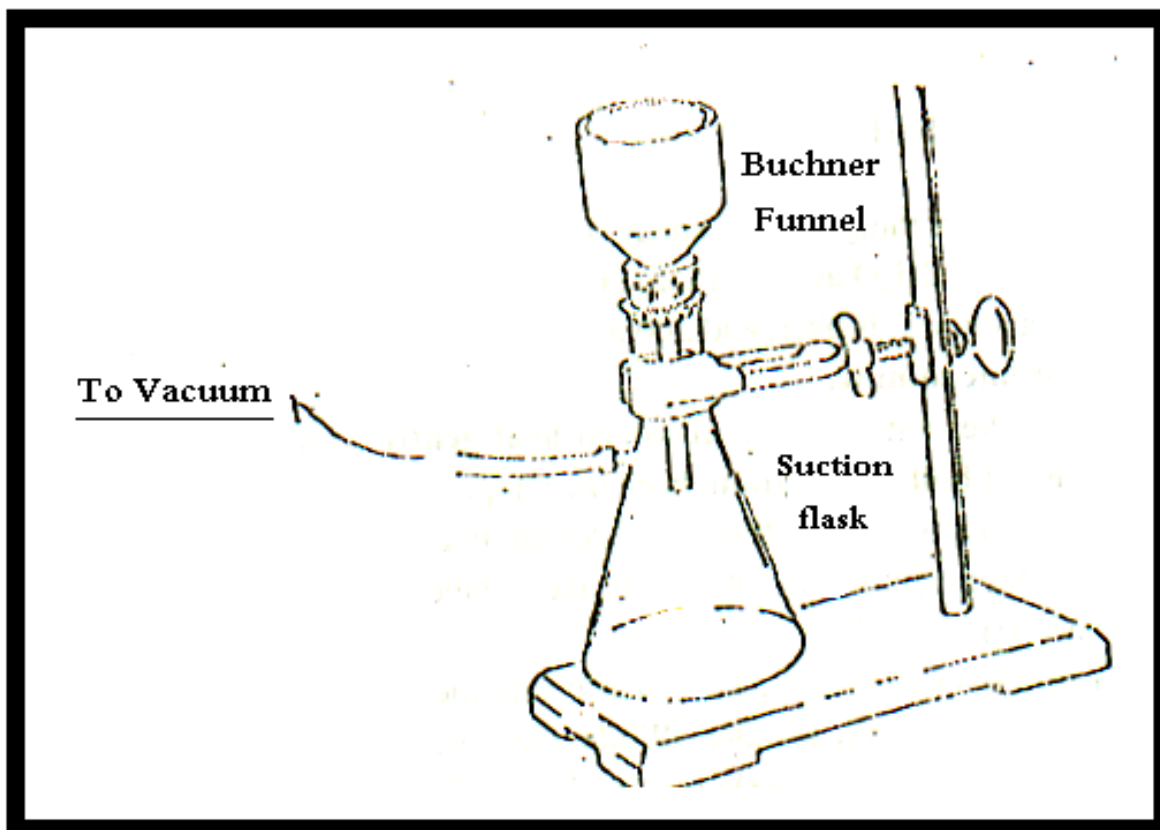


Fig. (1): Setup for vacuum filtration

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
 ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
 Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

### Experiment (3)

### Limiting Reactant

### Report Sheet

#### Part (1):- Preparation of $Ba_3(PO_4)_2$

Unknown number \_\_\_\_\_

<u>Data:</u>	<b>Trial 1</b>	<b>Trial 2</b>
1. Mass of mixture (g)	_____	_____
2. Mass of filter Paper (g)	_____	_____
3. Mass of Filter paper + $Ba_3(PO_4)_2$ (g)	_____	_____
4. Mass of $Ba_3(PO_4)_2$ (g)	_____	_____

#### Part (2): - Determination of the Limiting reactant

<u>Data:</u>	<b>Trial 1</b>	<b>Trial 2</b>
1. Limiting reactant in salt mixture	_____	_____
2. Excess reactant in salt mixture	_____	_____
<b>Calculations:</b>		
1. Mass of $Ba_3(PO_4)_2$ (g)	_____	_____
2. Moles $Ba_3(PO_4)_2$ : (FM = 602.2 g/mole)	_____	_____
3. Moles $BaCl_2 \cdot 2H_2O$ ; (FM = 244.2g/mole)	_____	_____
4. Mass of $BaCl_2 \cdot 2H_2O$ reacted (g)	_____	_____
5. Moles $Na_3PO_4 \cdot 12H_2O$ reacted; (FM = 380.2 g/mole)	_____	_____
6. Mass of $Na_3PO_4 \cdot 12H_2O$ reacted (g)	_____	_____
7. Mass of salt mixture (g)	_____	_____
8. Mass of excess reactant (g)	_____	_____
9. % Of $BaCl_2 \cdot 2H_2O$ in the mixture	_____	_____
10. % Of $Na_3PO_4 \cdot 12H_2O$ in the mixture	_____	_____

**Show your calculation on a separate paper.**



**Experiment (3)**  
**Limiting Reactant**

**PostLab**

1. What is the effect of heating the solution on the particle size of  $\text{Ba}_3(\text{PO}_4)_2$  precipitate?

---

---

---

2. Describe the effect of each of the following factors (whether increase, decrease, or has no effect) on the actual yield of  $\text{Ba}_3(\text{PO}_4)_2$

- a. Using a coarse filter paper instead of one with fine porosity.

---

---

- b. Insufficient washing of the precipitate.

---

- c. If the precipitate is washed with an acidic solution.

---

- d. If the precipitate was not dried completely.

---

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (4)

### Identification of a Compound: Physical Properties

#### Prelab

1. Define:-

A. Solubility:

---

---

---

B. Boiling point:

---

---

---

2. A student's liquid unknown boils at approximately 69°C, is insoluble in water but soluble in cyclohexane. Its density is 0.65g / ml. Which chemical in table (1) is the unknown?

---

---

---

3. What physical property, measurable in this experiment, distinguishes *cyclohexane* from *cyclohexene*?

---

---

---

4. Using the apparatus described in this experiment, when should the boiling point of a liquid be recorded?

---

---

---

## Experiment (4)

### Identification of a Compound: Physical Properties

#### **Objective:**

- To identify a compound based on its physical properties
- To learn how to properly assemble a boiling point apparatus

#### **Introduction:**

Gold is yellow, roses are red, salt is white and crystalline, lead is "heavy" and water is a colorless clear liquid. These are all physical properties of chemical substances. Physical properties are often used for the identification of a chemical substance. The more common physical properties include color, odor, density, solubility, melting point, boiling point, and whether the substance is a solid, liquid, or gas. These properties can be observed without any knowledge of the reactivity of the compound. Additional tests, tests that reveal more information about its chemical or physical properties, may be necessary, however, to confirm its purity.

The purpose of this experiment is to identify a chemical compound from the data collected on its solubility, density, and boiling point.

#### **Solubility:**

The solubility of a compound is the maximum mass (usually in grams) of the substance that dissolves in a fixed mass (usually 100g) of solvent at a given temperature. A chemical has different solubilities in different solvents, reflecting the differences in the molecular composition of the solute from that of the solvent. For instance, some substances such as table salt, are soluble in water but insoluble in gasoline. In this experiment we will examine the solubility of a substance in three different solvents: water, cyclohexane, (C<sub>6</sub>H<sub>12</sub>) and ethanol (C<sub>2</sub>H<sub>5</sub>OH). Qualitatively each substance is recorded as being soluble (complete dissolution), slightly soluble (partial dissolution), or insoluble.

#### **Density:**

The density is the mass per unit volume. Substances with a large density have a large mass in a small volume. While we commonly say that Lead is "heavy" we really mean that lead has a high density-it doesn't require much volume in order to have a large mass of lead.

**Boiling point:**

When a liquid is gradually heated, there is a temperature at which bubbles form spontaneously and continue to form until the entire volume of the liquid has been converted to a gas. This "constant" temperature is called the boiling point of the liquid. At the boiling point the temperature of the escaping vapor equals that of the liquid; this temperature depends upon the prevailing atmospheric pressure. The boiling point values listed in Table (1) are measured at normal atmospheric pressure (called one atmosphere of pressure).

The boiling point of a substance is characteristic of the magnitude of the forces acting between molecules, called *intermolecular forces*. The greater the magnitude of the intermolecular forces, the higher will be boiling point. Since the magnitude of the intermolecular forces is different for each compound, each has a characteristic boiling point.

**Apparatus and Chemicals:**

75 mm test tube (4), dropper, 150 ml-beaker, pipette, thermometers 110°C, rubber band, and capillary tube. Beaker 250 ml, Burner, wire gauze, O-ring. Stirring rod, Ethanol, Cyclohexane, Acetone, and about 15ml of an unknown liquid.

## Experimental Procedure:

**Note:** (Ask your instructor for an unknown listed in Table 1)

### A. Solubility:

#### 1. Solubility in Cyclohexane:

Into a 75-mm test tube add 3 drops of your liquid unknown to 1 ml of cyclohexane. Agitate the contents (Figure 1) Does each drop dissolve? Partially dissolve? Record your observations.

#### 2. Solubility in water and Ethanol:

Repeat the test with water and ethanol as solvents. Describe the solubility of your unknown as soluble, slightly soluble, or insoluble in each solvent.

**Table (1) physical properties of some common laboratory chemicals**

**Symbols used: *i* = insoluble, *sls* = slightly soluble, *s* = soluble**

Compound	Density(g/ml)	Boiling Point(C)	Solubility		
			H <sub>2</sub> O	C <sub>6</sub> H <sub>12</sub>	C <sub>2</sub> H <sub>5</sub> OH
Acetone	0.79	56	s	s	s
2-butanone	0.805	80	s	s	s
Cyclohexane	0.79	80.74	i	-	s
Cyclohexene	0.81	83	i	s	s
Ethanol	0.79	79	sls	s	-
Ethylacetate	0.90	77	s	s	s
Heptane	0.684	98	i	s	s
n-hexane	0.66	68	i	s	s
1-hexene	0.67	63	i	s	s
Isopropanol	0.79	83	s	s	s
Methanol	0.79	65	s	s	s
n-propanol	0.805	97	s	s	s
Water	1.00	100	-	i	s

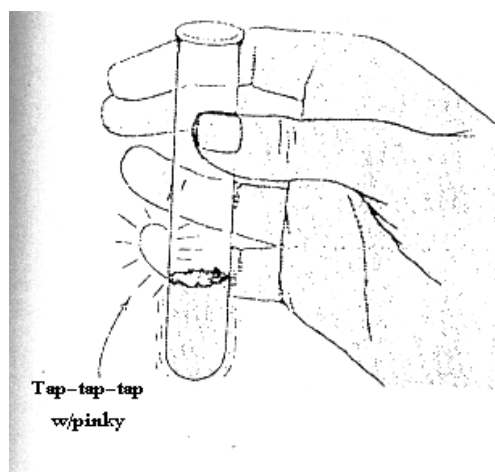


Figure (1) shake contents of the test tube

## B. Density:

Complete two trials to determine the density of your unknown.

Determine Mass and Volume of Unknown. Pipette 2ml of the liquid into a beaker and weigh. Calculate the density of the liquid.

## C. Boiling point:

Check out a 110°C thermometer. Assemble the apparatus in Figure 2. Complete two measurements for the boiling point of your unknown.

1. Assemble the Boiling Point Apparatus. Place 10 drops of the unknown liquid into a 75-mm test tube. Position it beside the thermometer bulb using a rubber band. Invert (open end down) a 10-cm capillary tube in the liquid. Place the apparatus into the water bath.
2. Measure the Boiling Point. Slowly heat the water in the water bath while stirring with the stirring rod. When a rapid and continuous stream of bubbles escapes the capillary tube, discontinue heating the water bath. The stream of bubbles slows as the bath cools. When the bubbles cease to escape and before the liquid re-enters the capillary tube, record the temperature. This is the boiling point of the liquid.
3. Repeat the Measurement. Determine the boiling point a second time. The same (volume of) liquid may be used. However, the capillary tube must be removed and re-inserted before heating resumes.

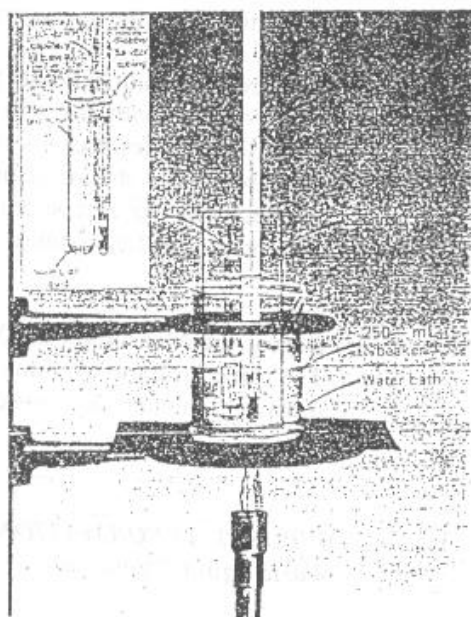


Figure (2) apparatus for determining the boiling point of a liquid

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_

ID Number: \_\_\_\_\_ Date : \_\_\_\_\_

Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

### **Experiment (4)**

## **Identification of a Compound: Physical Properties**

### **Report Sheet**

#### **Solubility:**

Unknown Number \_\_\_\_\_

	Soluble	Slightly Soluble	Insoluble
1. Solubility in cyclohexane	_____	_____	_____
2. Solubility in ethanol	_____	_____	_____
3. Solubility in water	_____	_____	_____

#### **Density:**

	Trial 1	Trial 2
1. Volume of liquid (ml)	_____	_____
2. Mass of liquid (g)	_____	_____
3. Density of liquid (g/ml)	_____	_____
4. Average density of liquid (g/m)	_____	

#### **Boiling point:**

1. Observed boiling point (c) \_\_\_\_\_

Name of unknown: \_\_\_\_\_

## Experiment (4)

### Identification of a Compound: Physical Properties

#### Post Lab

1. How does atmospheric pressure affect the boiling point of a liquid? Ask your instructor.

---

---

---

2. If several drops of liquid unknown (Part B) cling to the pipette wall after delivery will the density of the unknown be reported too high or too low? Explain

---

---

---

3. (A). In part C.2 if the boiling point is recorded when bubbles are rapidly escaping the capillary tube, will it be recorded too high or too low? Explain.

---

---

---

- (B). If the boiling point is recorded after the liquid enters the capillary tube (after the heat is removed), will it be recorded too high or too low? Explain

---

---

---



## Experiment (5)

### Tests for Cations and Anions

#### Objective:

- To observe reactions of several common cations and anions.
- To use chemical tests to identify an unknown salt.

#### Introduction:

One of the most common and important tasks for a chemist is to identify unknown compounds. In this experiment you will use simple chemical tests to identify an unknown salt containing one cation and one anion. We have limited the possible ions in your unknown to the following:

Cations	Anions
potassium, $K^+$	Sulfate, $SO_4^{2-}$
calcium, $Ca^{2+}$	Hydrogen carbonate, $HCO_3^-$
iron(III), $Fe^{3+}$	Chloride, $Cl^-$
ammonium, $NH_4^+$	Bromide, $Br^-$

Your unknown could be any of the 16 possible combinations of these cations and anions, for example  $KBr$ ,  $FeCl_3$ ,  $Ca(HCO_3)_2$ ,  $(NH_4)_2SO_4$ , etc.

#### Chemicals and equipment:

##### 0.1M solutions of:

sodium sulfate ,  $Na_2SO_4$   
potassium chloride ,  $KCl$   
sodium chloride ,  $NaCl$   
sodium bromide ,  $NaBr$   
sodium hydrogen carbonate,  $NaHCO_3$   
silver nitrate ,  $AgNO_3$   
barium chloride ,  $BaCl_2$   
sodium hydroxide ,  $NaOH$   
sodium oxalate ,  $Na_2C_2O_4$   
iron(III) sulfate ,  $Fe_2(SO_4)_3$   
potassium thiocyanate ,  $KSCN$   
calcium nitrate ,  $Ca(NO_3)_2$   
ammonium nitrate ,  $NH_4NO_3$   
2M nitric acid,  $HNO_3$   
2M hydrochloric acid,  $HCl$

Two test tube holders  
Bunsen burner  
Platinum wire at end of glass rod  
Red litmus paper  
Adhesive labels  
Unknown salt  
Test tube rack with test tubes

The procedure used in these tests is an example of **qualitative** analysis. You do not need to measure out volumes exactly using a measuring cylinder. In fact, doing so, would only waste time and increase the chance of contamination. Volumes of 1 or 2 mL are used. If you want to know how much this will fill a test tube, you can measure out 1 or 2 mL of distilled water in a measuring cylinder and pour it into a test tube. In the tests you can use approximately the same volumes of the solutions.

Use **clean** test tubes. Dirty glassware will give confusing or misleading results. After use, wash the test tubes several times with tap water and finally twice with distilled water from your wash bottle.

Each test is performed with a known compound and with your unknown salt in two test tubes held side by side so that the results can be compared. A positive test is a test that confirms the presence of a certain ion. When you perform the test on the known compound you will, of course, obtain a positive test. When you do the same test on your unknown you may see the same reaction, a different reaction or no reaction at all. Write down what you see using a complete sentence. For example: "A white precipitate formed", "The solution remains colorless", "A clear blue solution forms", etc. Do not write 'No change.' **You get a positive test for your unknown if you get the same result as was obtained for the known compound.** You should label each test tube to remember which tube contains the known compound and which contains the unknown.

To gain experience you should do all the tests for the cations and the anions even if you have already established the identity after a few tests. You can do the tests in any order.

## Experimental Procedure:

### A. Prepare a solution of your unknown

Place four spatula tips of your unknown sample in a 100 mL conical flask and add about 20 mL distilled water. Swirl to dissolve the salt. If you have difficulty in dissolving the salt, try to heat it in a hot water bath. Use this unknown solution in the following tests. More solution can be made up later if necessary.

### B. Tests for cations

#### 1. Test for calcium ion, $Ca^{2+}$

Take two clean test tubes. Place about 2 mL of calcium nitrate solution (known solution of calcium ion) in one test tube and about 2 mL of your unknown solution in the other test tube. Add about 10 drops of the sodium oxalate solution to each tube. Record your observation and conclusion below.

Solution	Observation	Conclusion
$Ca(NO_3)_2$		$Ca^{2+}$ is present
Unknown		

#### 2. Test for iron(III), $Fe^{3+}$

Take 2 mL of iron(III) sulfate in one test tube and 2 mL of your unknown solution in another test tube.

Add about 10 drops of potassium thiocyanate (KSCN) solution to each test tube. Shake to mix. Record your observations and conclusion below.

Solution	Observation	Conclusion
$Fe_2(SO_4)_3$		$Fe^{3+}$ is present
Unknown		

### 3. Test for ammonium ion, $\text{NH}_4^+$

Place 2 mL ammonium nitrate solution in one test tube and 2 mL unknown solution in another test tube.

Add 1 mL sodium hydroxide solution to each tube.

Warm each tube in a hot water bath for about one minute. Hold a moistened (wet) piece of red litmus paper at the mouth of each test tube, taking care not to touch the side of the tube as you are testing for any gas evolved.

Solution	Observation	Conclusion
$\text{NH}_4\text{NO}_3$		
Unknown		

#### Conclusion from the cation tests

Unknown No. \_\_\_\_\_

The cation \_\_\_\_\_

### C. Tests for anions

#### 1. Test for sulfate ion, $\text{SO}_4^{2-}$

In one test tube place 2 mL sodium sulfate solution and in another tube 2 mL of your unknown solution.

Add 10 drops of hydrochloric acid to each tube, shake gently and then add about 10 drops of barium chloride solution. Record your observations and conclusion below.

Solution	Observation	Conclusion
$\text{Na}_2\text{SO}_4$		
Unknown		

#### 2. Test for hydrogen carbonate ion, $\text{HCO}_3^-$

Place 2 mL sodium hydrogen carbonate solution in one test tube and 2 mL of your unknown in another test tube.

Add 2 mL hydrochloric acid to each tube and watch carefully for any **immediate** reaction. Record your observations and conclusion below.

Solution	Observation	Conclusion
$\text{NaHCO}_3$		
Unknown		

3. *Test for chloride and bromide ions,  $Cl^-$  and  $Br^-$*

Take **three** clean, labeled test tubes. Place 2 mL sodium chloride solution in one test tube, 2 mL sodium bromide solution in another test tube and 2 mL of your unknown in the third test tube.

Add about 10 drops of nitric acid and about 5 drops of silver nitrate solution to each tube. If you get a precipitate, look carefully at the color of the precipitate.

<b>Solution</b>	<b>Observation</b>	<b>Conclusion</b>
NaCl		
NaBr		
Unknown		

*Conclusion from the anion tests*

Unknown No. \_\_\_\_\_

The anion \_\_\_\_\_

**Final conclusion**

The formula of salt No. \_\_\_\_\_ is \_\_\_\_\_

The name of the compound is: \_\_\_\_\_

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## **Experiment (6)**

### **Molar Mass of a Volatile Liquid**

#### **Prelab**

- 1) For which of the following compounds can we determine its molar mass using the method described in this experiment? Give reasons. Benzene (b.p.78°C), Glycerol (b.p.180°C).

---

---

---

---

---

- 2) A cylinder contains compressed hydrogen gas and the mass of the hydrogen is 20.0 g. What mass of oxygen would be contained in an identical cylinder at the same temperature and pressure?

---

---

---

---

---

---

---

## Experiment (6)

### Molar Mass of a Volatile Liquid

**Objective:**

To determine the molar mass of an unknown volatile liquid.

**Introduction:**

According to Avogadro's principle, equal volumes of gases at the same temperature and pressure contain equal number of particles (molecules). Dumas used this principle in 1826 to determine atomic and molar masses. Dumas measured the densities of volatile liquids (easily vaporized), and from these densities he determined the atomic and molar masses of many gases. In this experiment the molar mass of a volatile liquid is determined using the modified Dumas apparatus and procedure. A volatile liquid is transferred to a weighed Erlenmeyer flask. The flask is then heated to the boiling point of water to vaporize the liquid and expel air and excess liquid. Finally the flask is cooled to room temperature. The flask with contents is weighed again and the mass of the condensed vapor is calculated. The molar mass is determined using the ideal gas law:

$$PV = nRT \quad (1)$$

The number of moles,  $n$ , is given by

$$n = (m/M) \quad (2)$$

The ideal gas equation becomes:

$$PV = (m/M) RT \quad (3)$$

Where:  $P$  = atmospheric pressure (atm)

$V$  = volume of vapor (L)

$m$  = mass of vapor (g)

$M$  = molar mass of volatile liquid (g/mol)

$R$  = gas constant =  $0.0821 \text{ atm L K}^{-1} \text{ mol}^{-1}$


$T$  = temperature of boiling water (K).

### Experimental Procedure:

1. Weigh together a clean and **dry** 125 ml Erlenmeyer flask, a copper wire and an aluminum foil (8 x 8 cm) on a Top-loading balance and record the weight as accurately as possible.
2. Using a graduated cylinder, pour about 10 ml sample of the volatile liquid (supplied by your instructor) into the Erlenmeyer flask.
3. Center the aluminum foil (8 x 8 cm) over the mouth of the flask and secure it in place with the copper wire.
4. Place a larger piece of aluminum foil (10 x 10 cm) over the first foil, cover and secure it in place with a rubber band.
5. Make a small hole through the aluminum foils with a sharp pin.
6. Place the Erlenmeyer flask into a 400 ml beaker and clamp it to a stand.
7. Place some water in the beaker, to cover up the liquid level in the flask, along with few boiling stones.
8. Place the beaker, with the flask over Bunsen burner as shown in the following figure.
9. Heat the beaker slowly till the water boils (this step should take about 10 minutes).
10. Allow the water to boil gently, till all the volatile liquid has vaporized. Continue gentle heating of water for another 5 minutes.
11. Measure the temperature of the boiling water and record it.
12. Turn off the Bunsen burner.
13. Using a towel, immediately remove the **Erlenmeyer** flask from the beaker.
14. Dry the outside of the flask and the aluminum foil **completely**.
15. Remove the outer aluminum foil **and rubber band** only.
16. Allow the Erlenmeyer flask to sit for 15 minutes to cool to room temperature.
17. While waiting for the Erlenmeyer flask to cool, read the atmospheric pressure on the laboratory barometer.
18. Weigh the Erlenmeyer flask, with aluminum cover and wire still on, accurately and record the mass.



19. Determine the volume of the flask by filling it to the rim with water and measuring the volume of the water with a graduated cylinder. Record the volume.
20. Repeat the above steps for a second trial.

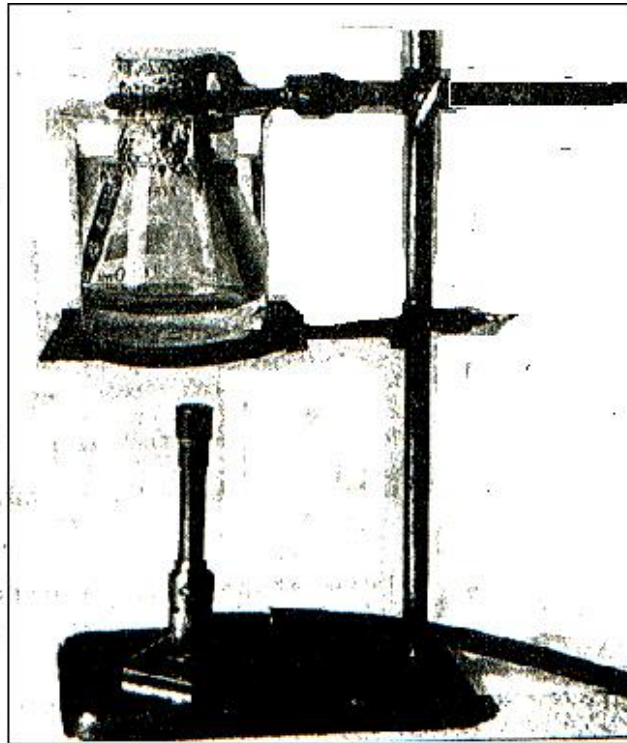

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (6)

### Molar Mass of a Volatile Liquid

#### Report Sheet

Unknown No. \_\_\_\_\_

#### Data:

	<u>Trial I</u>	<u>Trial II</u>
1. Mass of empty Erlenmeyer flask + Aluminum foil + Wire (g)	_____	_____
2. Temperature of boiling water (K)	_____	_____
3. Atmospheric pressure (atm)	_____	_____
4. Mass of Erlenmeyer flask + Condensed vapor + aluminum foil + wire (g)	_____	_____
5. Mass of condensed vapor (g)	_____	_____
6. Volume of Erlenmeyer flask	_____	_____
7. Molar mass of the liquid, calculated	_____	_____
8. Average molar mass	_____	

**Note:** *A more accurate calculation of the molar mass uses a corrected mass for the vapor which takes into consideration the mass of air displaced by vapor at room temperature.*

## Experiment (6)

### Molar Mass of a Volatile Liquid

#### PostLab

1. A flask weighs 40.1305g when clean, dry, evacuated, 138.2410g when filled with water (density = 0.9970 g/ml) and 40.2487g when filled with a gaseous substance at 470.4 torr and 96°C. What is the molar mass (g/mol) of the gas? [R = 0.082 L.atm/mol. K]

---

---

---

---

---

---

---

---

---

---

---

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## **Experiment (7)**

### **Determination of the Molar Volume of Hydrogen Gas**

#### **PreLab**

A student at The Hashemite University wants to determine experimentally the volume occupied by one mole of  $H_2$  gas at STP. She reacts 0.1471 g of Zn with excess  $HCl(aq)$  and collects 56.09 ml of gas over water at 22 °C and 757.8 torr. The vapor pressure of water at 22 °C is 19.8 torr.

1. Use the data given above to calculate:
  - (i) The volume occupied by one mole of dry  $H_2$  at 22 °C and 760 torr.
  - (ii) The volume occupied by one mole of dry  $H_2$  at STP.

Note: do NOT simply use the ideal gas law equation!

2. Name the gas laws which you used in your calculations.

## Experiment (7)

### Determination of the Molar Volume of Hydrogen Gas

#### Objectives:

- To determine the volume, in liters, of one mole of dry hydrogen ( $H_2$ ) gas at STP

#### Introduction:

Hydrogen gas, or dihydrogen, ( $H_2$ ) can be produced from the reaction of active metals with acids. In this experiment you will generate  $H_2$  by reacting magnesium metal with hydrochloric acid. You will collect the  $H_2$  in an inverted burette by water displacement and measure the volume of the gas. You can then calculate the volume of dry gas at STP if you know the room temperature, the barometric pressure, and the vapor pressure of water. We know that the ideal gas law predicts a volume of  $22.4 \text{ L mol}^{-1}$ , but in this experiment we want to *measure* this for ourselves.

#### Apparatus and Chemicals:

50 mL burette, one- or two-hole stopper to fit burette, ring stand with clamp, 600 mL beaker, 10 mL and 1000 mL graduated cylinders, fine copper wire, 6 M HCl, 5-6 cm piece of magnesium ribbon.

#### Experimental Procedure:

1. The first step is to record the volume contained in your burette between the 50 mL mark and the sealed end. It is written on your burette
2. Weigh a piece of magnesium ribbon with a mass of approximately 0.035–0.04 g on an analytical balance. Record its mass to four decimal places.
3. Fold the magnesium ribbon into a small pellet. Wrap it tightly all around with copper wire so that no small pieces can escape as it dissolves. Leave about 5 cm of copper wire to serve as a handle. See Figure 1.
4. Clamp your burette in an upright position in the ring stand. Fill a 600 mL beaker about two-thirds full of tap water and place near the ring stand.
5. Measure about 10 mL of 6 M HCl in a graduated cylinder. Tilt the burette at an angle from the upright position and pour in the HCl.
6. With the burette in the same position, SLOWLY fill it with tap water from the beaker. While pouring, rinse any acid that may be on the sides of the tube so that the liquid in the top of the tube will contain very little acid. TRY TO AVOID STIRRING UP THE ACID LAYER IN THE BOTTOM OF THE TUBE.
7. While holding the copper coil by the handle, insert the metal about 3 cm down into the burette. Hook the copper wire over the edge of the tube and hold it there by inserting the stopper. Fill the burette completely so that the stopper displaces a little water when put in place. See the left of Figure 1.

8. Cover the hole in the stopper with your finger and invert the burette in the 600 mL beaker of water, as shown in the middle of Figure 1. Clamp the burette in place. The acid being more dense than water will flow down to the bottom and eventually react with the metal.
9. After the reaction stops, wait for about 5 minutes to allow the burette to come to room temperature. If any bubbles stick to the sides of the burette, tap them gently.
10. Cover the hole in the stopper with your finger and transfer the burette to a large graduated cylinder almost full of room temperature water. See Figure 2. Raise or lower the burette until the level of the liquid inside the burette is the same as the level outside the burette. This permits you to measure the volume of the gases in the burette ( $H_2$  and water vapor) at room pressure.
11. Read the burette. The calibrations are upside down, so make sure you read it correctly. Ask for help if you are not sure of the reading. Record this reading on your data sheet.
12. Remove the burette from the water and pour the acid solution it contains down the sink. Rinse the burette with tap water and then with distilled water.
13. Record the temperature of the water in the large graduated cylinder, the room temperature, and the room pressure as given by your instructor.
14. Do the experiment a second time.

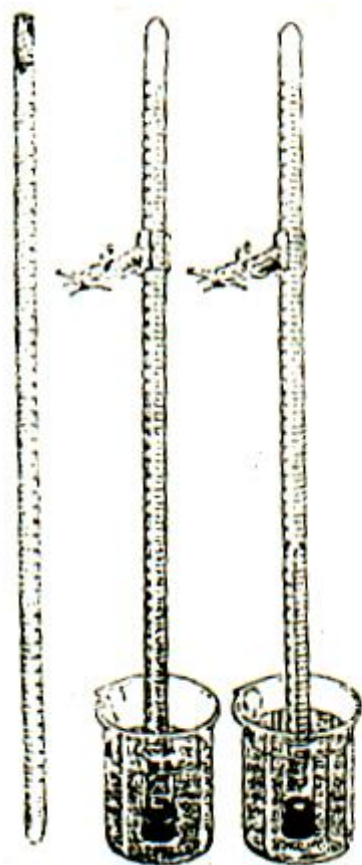


Figure 1 Manipulating the burette.



Figure 2 Measuring the volume of gas.

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
 ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
 Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (7)

### Determination of the Molar Volume of Hydrogen Gas

#### Report Sheet

	Trial 1	Trial 2
Volume of burette from 50 mL mark to the sealed end		
Mass of piece of magnesium ribbon		
Burette reading from Step 11		
Temperature of the water in the graduated cylinder (Step 13)		
Room temperature		
Room pressure		
Vapor pressure of water—use the water temperature from Step 13 and the Table below.		

Vapor Pressure of Water at Various Temperatures			
Temperature °C	Pressure torr	Temperature °C	Pressure torr
16	13.6	21	18.6
17	14.5	22	19.8
18	15.5	23	21.0
19	16.5	24	22.4
20	17.5	25	23.8

### Calculations:

Set up each problem completely and **SHOW UNITS**. Cancel the units when you do the calculation. Use the correct number of significant figures in your answers.

1. Write the balanced equation for the reaction.

---

2. Determine the number of moles of magnesium used. Find the average for the 2 trials.

---

---

3. Determine the volume of gas produced. Find this by the following equation:

Volume of gas = (Volume of burette from 50 mL to sealed end) + (50.0 mL – burette reading). Find the average for the 2 trials.

*Use the average of your data from the previous page for these calculations.*

---

---

4. Determine the partial pressure of the H<sub>2</sub> gas. Since the H<sub>2</sub> is collected over water, the gas in the burette consists of a mixture of H<sub>2</sub> and water vapor. The total pressure is  $P_{\text{room}} = P_{\text{H}_2} + P_{\text{H}_2\text{O}}$ .

---

---

---

5. Using the results of your calculations in the previous steps, determine the volume (in liters) of the H<sub>2</sub> at one atmosphere pressure (760 torr) if no water were present.

---

---

---

6. Using the results of your calculations in the previous steps, calculate the molar volume (in liters/mol) of dry H<sub>2</sub> at room temperature and 1 atmosphere pressure. (*Don't* simply use the ideal gas law, use your own results!)

---

---

7. Using the results of your calculations in the previous step, calculate the molar volume (in liters/mol) of dry H<sub>2</sub> at STP.



Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## **Experiment (8)**

### **Colligative Properties: Molar Mass Determination**

#### **PreLab**

1. Students prepared two cyclohexane solutions having the same mass of solute. However Student 1 used 13 g of cyclohexane, Student 2 used 15 g. Which student will observe the larger freezing point change? Explain.

---

---

---

2. A 0.597 g sample of a non-electrolyte dissolves in 20.0g of cyclohexane. The freezing point depression is 3.62 °C. What is the molar mass of the non-electrolyte? ( $K_f$  for Cyclohexane is 20.0 °C Kg/mol).

---

---

---

---

---

---

---

---

## Experiment (8)

### Colligative Properties: Molar Mass Determination

#### Objective:

To determine the molar mass of a non-volatile, non-electrolyte by observing the difference between the freezing points of a solvent and a solution.

#### Introduction:

The addition of a non-volatile solute to a solvent produces characteristic changes in a solvent's physical properties. These changes are: **vapor pressure lowering**, **boiling point elevation** and **freezing point depression**. For example when antifreeze (ethylene glycol) is added to an automobile's cooling system it will prevent freeze-up in the winter and boiling in the summer because the antifreeze-water solution has a lower freezing point and a higher boiling point than pure water.

Vapor pressure lowering, boiling point elevation, and freezing point depression are called **colligative properties** of solutions. These properties are governed by the number, rather than the type, of solute particles dissolved in the solvent.

The change in colligative properties of a solution is directly proportional to the amount of solute dissolved in the solvent. The freezing point depression,  $\Delta T_f$ , and the boiling point elevation,  $\Delta T_b$ , are proportional to the molality,  $m$ , of the solute in the solution:

$$\Delta T_f = K_f m \quad (1)$$

$$\Delta T_b = K_b m \quad (2)$$

Where

$K_f$  and  $K_b$  are the molar freezing point and boiling point constants for the solvent.

$$m = \frac{\text{mol solute}}{\text{Kg solvent}} = \frac{(\text{g / mol. wt.}) \text{ solute}}{\text{Kg solvent}} \quad (3)$$

In this experiment you will determine the molar mass of a non-volatile, non-electrolyte (solute) in cyclohexane by measuring the freezing points for the pure cyclohexane and the solution, as well as the mass of solute and solvent. From  $K_f$  and equations (1) and (3), you will calculate the formula weight of the solute.

Cyclohexane's freezing point and that for the solutions are obtained from a cooling curve; a graph representing a decreasing temperature as a function of time (Figure 1). Cyclohexane's cooling curve reaches a plateau at its freezing point; extrapolation of the plateau to the temperature axis determines its freezing point. The solution's cooling curve does not reach a plateau, but continues to decrease slowly as the cyclohexane freezes. Its freezing point is determined at the intersection of two lines drawn tangent to the curves above and below the freezing point (See Figure 1).

## Experimental Procedure:

### A. Freezing Point of Solvent:

1. Assemble the apparatus shown in Figure (2). Weigh ( $\pm 0.01$  g) a clean, dry 200 mm test tube. Add approximately 10 g of cyclohexane to the test tube and reweigh ( $\pm 0.01$  g) the tube and contents. Prepare about 300 mL of an ice-water slurry in a 400-mL beaker. Place the test tube and cyclohexane in the ice-water bath (Figure 2). Insert a thermometer, through one hole of a two-hole stopper, into the solvent to measure the temperature.
2. While stirring with the wire stirrer inserted through the stopper's other hole, record temperature ( $\pm 0.1^\circ\text{C}$ ) readings at timed intervals (30 sec) on the Report Sheet. The temperature remains virtually constant at the freezing point until solidification is almost complete. Continue recording until the temperature begins to drop again. On a linear graph paper (provided by your laboratory technician), plot the temperature ( $^\circ\text{C}$ , vertical axis) versus time (t, horizontal axis) to obtain cyclohexane's "cooling curve". The temperature corresponding to the plateau of the cooling curve is cyclohexane's freezing point.

### B. Freezing Point of Solvent-Unknown Solute:

Two freezing point trials for the solution are to be completed. The same unknown solute-cyclohexane solution can be used to repeat Parts B.1 and B.2.

1. Dry the outside of the test tube containing the cyclohexane and reweigh ( $\pm 0.01$  g). Ask the instructor for a solid unknown sample. Record your unknown number on your data sheet. Weigh ( $\pm 0.01$  g) approximately 0.2g of the unknown on weighing boat. Transfer the unknown to the cyclohexane in the 200-mm test tube. In the transfer be certain that **none** of the solid solute remains in the boat or adheres to the test tube's side. If some does, roll the test tube until the solute dissolves. If an unknown is a liquid, weigh ( $\pm 0.01$  g) about 3 mL of it in a 10-mL graduated cylinder, pour this into the cyclohexane; then reweigh the empty cylinder. The difference is the weight of the unknown liquid added to the cyclohexane.
2. Determine the solution's freezing point in the same way as the solvent's was found in Part A. When the solution is near the freezing point, record the temperature ( $\pm 0.1^\circ\text{C}$ ) at more frequent intervals (20 sec). Plot the temperature vs. time (on the same graph as in Part A) to get the solution's cooling curve. It will show a "break" at the temperature where freezing began. The break may not be as sharp as that for the pure solvent. Draw lines tangent to the curves above and below the freezing point; the lines' intersection is the solution's freezing point. Repeat the freezing point determination one more time (warm the tube and the solution to a temperature higher than the freezing point of the solution then cool again as before.) Have your instructor approve your graphs. After completing the collection of the experimental data, discard the solution into a bottle marked "waste solvent."

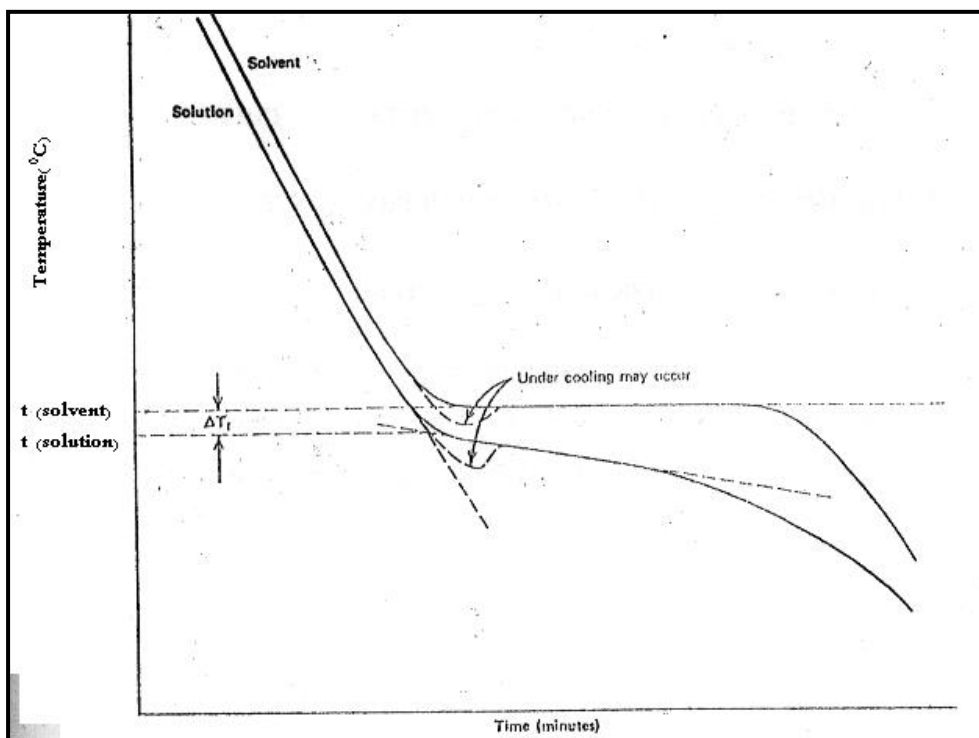


FIGURE 1. Cooling curves for solvent and solution

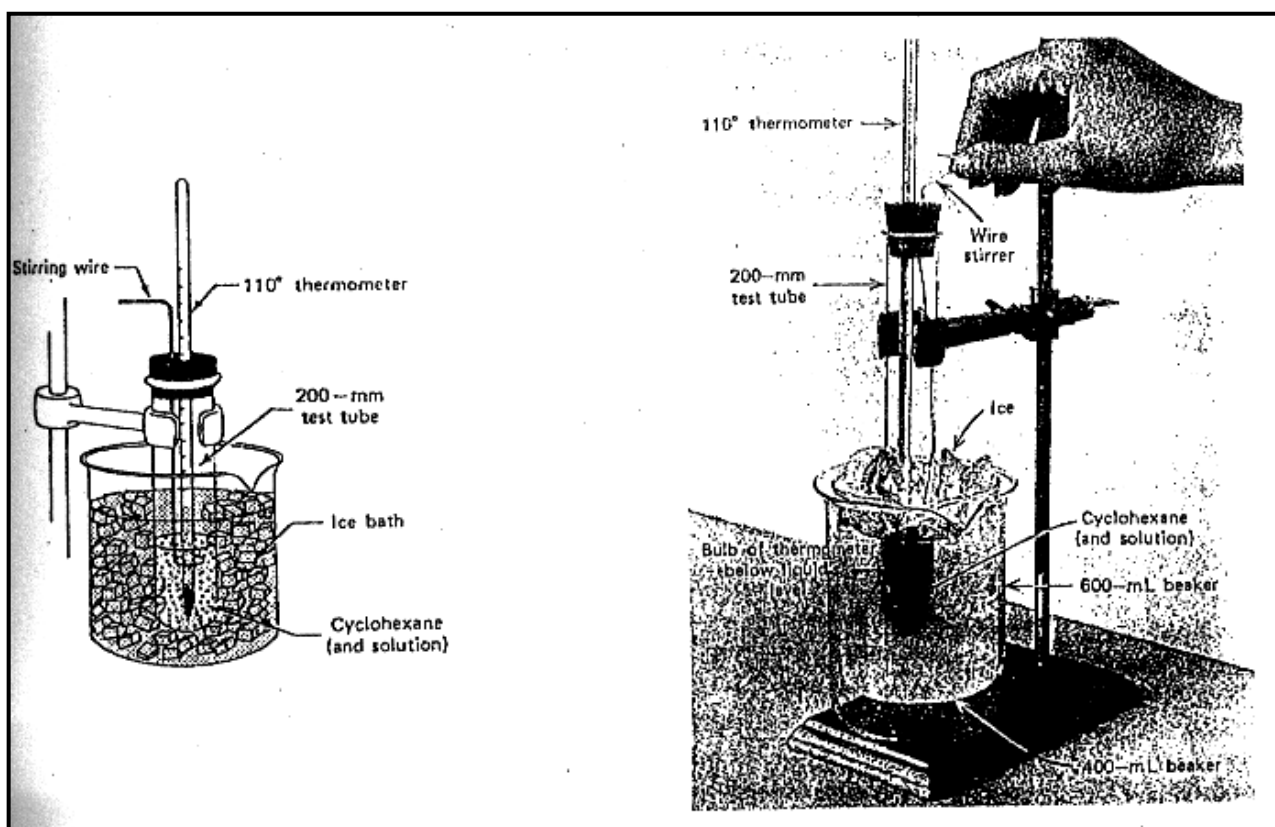


FIGURE 2. Freezing point apparatus

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (8)

### Colligative Properties: Molar Mass Determination

#### Report Sheet

Unknown Number \_\_\_\_\_

#### Data:

##### A. Freezing Point of Cyclohexane

1. Mass of test tube (g) \_\_\_\_\_
2. Mass of cyclohexane (g) \_\_\_\_\_
3. Freezing point, from cooling curve ( $^{\circ}\text{C}$ ) \_\_\_\_\_

##### B. Freezing Point of a Solvent-Unknown Solute

**Trial I**      **Trial II**

- |  |       |       |
|--|-------|-------|
| 1. Mass of test tube + cyclohexane (g)                       | _____ | _____ |
| 2. Mass of cyclohexane (g)                                   | _____ | _____ |
| 3. Mass of solute (g)  | _____ | _____ |
| 4. Freezing point, from cooling curve ( $^{\circ}\text{C}$ ) | _____ | _____ |

#### Calculations:

- |   |                                |
|---|--------------------------------|
| 1. Freezing point constant, $K_f$ , of cyclohexane            | 20.0 $^{\circ}\text{C.kg/mol}$ |
| 2. Freezing point change, $\Delta T_f$ ( $^{\circ}\text{C}$ ) | _____                          |
| 3. Mass of solute in solution (g)                             | _____                          |
| 4. Mass of cyclohexane in solution (kg)                       | _____                          |
| 5. Molar mass (g/mol)   | _____                          |
| 6. Average molar mass (g/mol)                                 | _____                          |

<b>DATA FOR COOLING CURVES</b>					
<b>PART A (SOLVENT)</b>		<b>PART B (SOLUTION)</b>			
<b>Trial 1</b>		<b>Trial 1</b>		<b>Trial 2</b>	
<b>Time(second</b>	<b>Temp.</b>	<b>Time</b>	<b>Temp.</b>	<b>Time</b>	<b>Temp.</b>
<b>0</b>		<b>0</b>		<b>0</b>	
<b>30</b>		<b>20</b>		<b>20</b>	
<b>60</b>		<b>40</b>		<b>40</b>	
<b>90</b>		<b>60</b>		<b>60</b>	
<b>120</b>		<b>80</b>		<b>80</b>	
<b>150</b>		<b>100</b>		<b>100</b>	
<b>180</b>		<b>120</b>		<b>120</b>	
<b>210</b>		<b>140</b>		<b>140</b>	
<b>240</b>		<b>160</b>		<b>160</b>	
<b>270</b>		<b>180</b>		<b>180</b>	
<b>300</b>		<b>200</b>		<b>200</b>	
<b>330</b>		<b>220</b>		<b>220</b>	
<b>360</b>		<b>240</b>		<b>240</b>	
<b>390</b>		<b>260</b>		<b>260</b>	
<b>420</b>		<b>280</b>		<b>280</b>	
<b>450</b>		<b>300</b>		<b>300</b>	
<b>480</b>		<b>320</b>		<b>320</b>	
<b>510</b>		<b>340</b>		<b>340</b>	
<b>540</b>		<b>360</b>		<b>360</b>	
<b>570</b>		<b>380</b>		<b>380</b>	
<b>600</b>		<b>400</b>		<b>400</b>	
<b>630</b>		<b>420</b>		<b>420</b>	
<b>660</b>		<b>440</b>		<b>440</b>	

- If needed, continue recording your data on separate sheets of papers and submit them with your Lab Report.

## Experiment (8)

### Colligative Properties: Molar Mass Determination

#### PostLab

1. If the solution's freezing point is erroneously read  $0.2^{\circ}\text{C}$  lower than it should be, will the unknown's calculated molar mass be too high or too low ? Explain!

---

---


2. How will the freezing point change of cyclohexane be affected ( compared with the freezing point change by a non-volatile, non-associating and non-dissociating solute) by:

- (a) A non-volatile solute that dissociates? Explain.

---

---

- (b) Two solutes that react according to the equation,

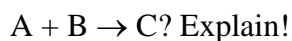

---

---

3. If some solute adheres to the test tube's wall in Part B.1, is the freezing point change greater or less than it should be? Explain!

---

---

---



## Experiment (9)

# Calorimetry


### Objective :

To measure heat of reaction

### Introduction :


This experiment is designed to measure the HEATS OF REACTION of a series of chemical reactions. There are several kinds of "heats of reactions", such as heat of solution, heat of neutralization, heat of combustion, heat of formation, etc. However, we shall study only the first two of these.

1. **Heat of solution**: it is the heat flow during the process of solution; as an example when NaOH dissolves in water, hydrated ions are formed and heat is released. The solution absorbs this heat and, therefore increases in temperature.

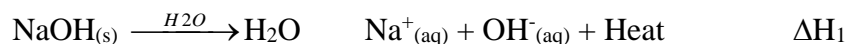

Where the heat of solution in kJ is  $\Delta H_1$

2. **Combined heat of solution and heat of neutralization**: in a second reaction, solid sodium hydroxide is reacted with an aqueous solution of hydrogen chloride. During this reaction, water and a dilute aqueous solution of sodium ions and chloride ions are formed. This time both the heat of solution and the heat of neutralization are released:

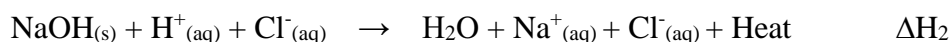

Where the combined heat of solution and heat of neutralization in kJ is  $\Delta H_2$ .

3. **Heat of neutralization**: It is the heat evolved, at constant pressure, when acids are neutralized with bases. **It** can be obtained from the two steps above. The difference between the **heat of reaction** obtained by dissolving solid NaOH in aqueous HCl solution (step 2) and the **heat of solution** obtained by dissolving solid NaOH in water (step 1).
4. The final experiment of this laboratory class will be to determine the heat of solution of an unknown salt, as in (step 1) above.

### Example 1

A 2.00 g sample of solid CsOH is dissolved in 200.0 mL of water in a calorimeter. The temperature of the water was raised from 22.3 to 23.4 °C. Calculate the heat of solution in kJ/mol. (assume the specific heat of the solution to be 4.184 J/g °C and the density of the solution to be 1.00 g/mL).

- Temperature change,  $\Delta T = 23.4 - 22.3 = 1.1\text{ }^\circ\text{C}$
- Total mass of solution, m, (200.0 mL water) (1.00 g/mL) + 2.00 g CsOH = 202.0 g
- Heat absorbed by the solution = heat evolved by CsOH =  
Mass of solution x specific heat of solution x  $\Delta T$  =  
 $202.0\text{ g} \times 4.184\text{ J/g }^\circ\text{C} \times 1.1^\circ\text{C} = 930\text{ J}$
- Moles of CsOH = mass/Molar mass =  $2.00\text{ g} / 149.9\text{ g/mol} = 0.01334\text{ mol}$
- Heat evolved per mole of CsOH =  $930\text{ J} / 0.01334\text{ mol} = -7.0 \times 10^4\text{ J/mol}$
- Heat of solution in KJ/mole =  $(-7.0 \times 10^4\text{ J/mol}) \times (1\text{ kJ}/10^3\text{ J}) = -70\text{ kJ/mol}$

Note: the minus sign shows that the reaction is exothermic; i.e heat is given off (evolved) by the system.

### Example 2

A 2.00 g sample of solid CsOH reacted with 200.0 mL of an aqueous solution of hydrogen chloride (HCl) in a calorimeter, the temperature of the solution increased from 22.3 to 24.3 °C. Calculate the heat of reaction in kJ/mol? (Assume the specific heat of the solution to be 4.184 J/g °C and the density of the solution to be 1.00 g/mL)

- Temperature change,  $\Delta T = 2.0\text{ }^\circ\text{C}$
- Total mass of solution = 202.0g
- Heat absorbed by the solution =  $202.0 \times 4.184 \times 2.0$   
 $= 1.7 \times 10^3\text{ J}$
- Moles of CsOH = 0.01334 mol
- Heat evolved per mole of CsOH =  $1.3 \times 10^5\text{ J/mol}$
- Heat of reaction in kJ/mol =  $1.3 \times 10^2\text{ kJ/mol}$

Note: The difference in the heat of reaction (Example 2) and the heat of solution (Example 1) is the **heat of neutralization**. So the heat of neutralization of CsOH with HCl is 60 KJ/mol

## **Experimental Procedure:**

Assume that the **density** of all solutions is 1.00 g/mL and their **specific heat** is 4.184 J/g °C .

### **A. Determination of the heat of solution for NaOH ( $\Delta H_1$ )**

1. Add 100.0 mL of water to the calorimeter.
2. Close the calorimeter and stir the water until a constant temperature is reached.  
Record the temperature of water.
3. Weigh accurately about 2.0 g of NaOH. You should weigh the NaOH fast, since it is very HYGROSCOPIC.
4. Add the solid NaOH to the water in the calorimeter. Stir gently. Read the temperature every 30 seconds. Record the maximum temperature reached.

### **B. Determination of the heat of reaction ( $\Delta H_2$ )**

1. Add 100.0 mL of 0.500 M HCl solution to the calorimeter.
2. Stir the solution as in procedure 1 above, and record the initial temperature.
3. Weigh accurately about 2.0 g of NaOH. You should weigh the NaOH fast, since it is very HYGROSCOPIC.
4. Add the solid NaOH to the HCl solution in the calorimeter and stir gently until the maximum temperature is reached. Record the maximum (final) temperature.

### **C. Determination of the heat of solution of an unknown salt:**

1. Pick up an unknown salt from the instructor and record the unknown number on the report sheet. You will be given the molar mass of the salt.
2. Repeat the procedure in part A, using 3.0 g of the unknown salt in place of NaOH.

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

**Experiment (9)**  
**Calorimetry**  
**Report Sheet**

**Part A: Determination of the heat of solution of NaOH ( $\Delta H_1$ )**

**Data:**

	<b>Trial I</b>	<b>Trial II</b>
1. Volume of water	_____	_____
2. Mass of water	_____	_____
3. Mass of NaOH	_____	_____
4. Molar mass of NaOH	40.0g/mol	
5. Moles of NaOH	_____	_____
6. Mass of solution	_____	_____
7. Initial temperature	_____	_____
8. Final temperature	_____	_____
9. $\Delta H_1$ of solution	_____	_____
10. Average $\Delta H_1$	_____	_____

**Part B: Determination of the heat of reaction of NaOH with HCl ( $\Delta H_2$ ) and the heat of neutralization ( $\Delta H_n$ )**

**Data:**

	<b>Trial I</b>	<b>Trial II</b>
1. Volume of HCl solution	_____	_____
2. Mass of HCl solution	_____	_____
3. Mass of NaOH	_____	_____
4. Molar mass of NaOH	40.0 g/mol	
5. Moles of NaOH reacted	_____	_____
6. Mass of solution	_____	_____
7. Initial temperature	_____	_____
8. Final temperature	_____	_____
9. $\Delta H_2$ of solution	_____	_____
10. Heat of neutralization ( $\Delta H_2 - \Delta H_1$ )	_____	_____
11. Average heat of neutralization	_____	

**Part C: Determination of the heat of solution of an unknown salt:**

**Data:**

	<b>Trial I</b>	<b>Trial II</b>
1. Volume of water	_____	_____
2. Mass of water	_____	_____
3. Mass of unknown salt	_____	_____
4. Initial temperature	_____	_____
5. Final temperature	_____	_____
6. $\Delta H$ of solution	_____	_____
7. Average $\Delta H$	_____	

**Show your calculation for Trial I:**

## Experiment (10)

### Electrochemistry


#### Objectives:

- To measure the cell potential for various electrochemical cells.
- To study the dependence of the cell potential concentration on the cell.



#### Introduction:

##### A. Electrochemical Cells:

An oxidation-reduction reaction (redox reaction) such as:

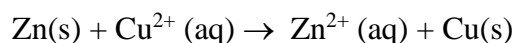

can be separated into two half-reactions, called half-cell reactions.

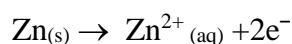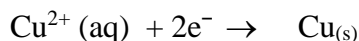

The first half-reaction is the reduction half-cell and the second half-reaction is the oxidation half-cell.


The above redox reaction can be carried out in a way such that the components of the oxidation and reduction half-cells are not mixed together, and the potential of the reaction can be measured, as shown in the Figure below.

Such a cell can be represented by:




The measured potential for this cell is 1.10 V. This value is referred to as the standard cell potential,  $E^{\circ}$  which indicates the voltage observed when all species are in their standard states (1.0 M for species in solution, and 1.0 atm. pressure for gases). Most commonly  $E^{\circ}$  values are measured at 298 K.

This cell potential,  $E^\circ$ , can be calculated from the two half-cell potentials:

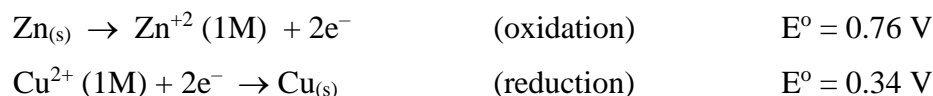

Net reaction:

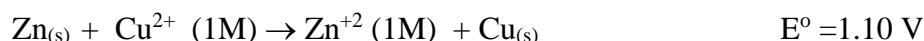

When a cell reaction is written so that  $E^\circ$  is positive, the reaction will proceed spontaneously from left to right when all components are in their standard states.

## **B. Concentration Cells**


In this experiment you will observe the effect of varying the concentration of a species, participating in an electrochemical oxidation-reduction reaction, on the cell potential. Since the potential of an electrode depends on the concentration of the ions used in the electrode half-cell, a cell may be constructed from two half-cells differing only in the concentration of ions. The principle behind this is simple, dilution of a concentrated solution is spontaneous, so this process gives a cell with positive  $E$ . For the process:



a cell can be constructed such as:



which represents the reaction:

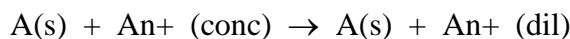

Since the two half-reactions are the same,  $E_o(\text{anode}) = E_o(\text{cath})$  so that  $E_o (\text{cell}) = 0$ .

Applying Nernst equation, at  $25^\circ\text{C}$ , we get:

$$E = -\frac{0.0592}{n} \log \frac{[\text{A}^{n+}]_{\text{dil}}}{[\text{A}^{n+}]_{\text{conc}}}$$

where  $n$  is the number of electrons transferred.

## **Experimental Procedure:**

### **A. Electrochemical Cells:**

1. Make a complete cell as shown in the Figure. Place 30 mL of 0.10 M  $\text{Zn}(\text{NO}_3)_2$  solution and 30 mL of 0.10 M  $\text{Cu}(\text{NO}_3)_2$  solution in separate 100 mL (or 150 mL) beakers.

Obtain a zinc plate and a copper plate from the store room, clean both sides of each plate (electrode) with a strip of sand paper.

2. Fill a U tube with saturated  $\text{KNO}_3$  solution, insert loose cotton plugs soaked with a saturated  $\text{KNO}_3$  solution in the ends of the tube leaving *no air bubbles*. Place the U tube across the beakers as a salt bridge.
3. Insert the metal plates in the solutions of their own ions.
4. Attach the appropriate electrodes to the voltmeter's terminals, using connecting wires, and read the voltage.
5. Remove the wires connecting the cell to the voltmeter and short-circuit the cell with a single piece of wire.

What chemical reaction is taking place?

Which way are the electrons flowing?

What changes can you observe which confirm that a chemical reaction is taking place?

6. Repeat steps 1-5 above making the following cells: Pb/Cu, and Pb/Zn. Record the voltage reading for each cell. (The standard reduction potential of  $\text{Pb}^{+2} + 2\text{e}^- \rightarrow \text{Pb}$  is  $-0.13 \text{ V}$ )

### **B. Concentration Cells**

Set up a simple concentration cell by using a Cu/0.10 M  $\text{Cu}(\text{NO}_3)_2$  half-cell as that prepared in section A, and a second half-cell prepared as follows:

1. Take 30 mL (the same volume as before) of 0.1M  $\text{KNO}_3$  solution and add to it exactly one drop of 0.10M  $\text{Cu}(\text{NO}_3)_2$  solution.
2. Clean a Cu strip with a sand paper and insert it into the beaker.
3. Connect the two half-cells by a salt bridge and measure the cell voltage.
4. Now add 10 drops of 0.10M  $\text{Na}_2\text{S}$  solution to the second, more dilute, half-cell. Since CuS is very insoluble the addition of  $\text{Na}_2\text{S}$  has the effect of decreasing  $\text{Cu}^{2+}$  concentration to an extremely small value.
5. Measure the cell voltage again.
6. Label the electrode signs and processes (draw a diagram).
7. Measure the temperature of the solution.



Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (10)

### Electrochemistry

### Report Sheet

#### Data and Calculations

##### A. Electrochemical Cells:

For each cell you made calculate  $E_o$  from your textbook table of electrode potentials. Compare your calculated values with those observed. Explain differences if any.

<u>Cell</u>	<u>Measured</u>	<u>Calculated</u>
Zn/Cu	_____	_____
Pb/Cu	_____	_____
Zn/Pb	_____	_____

Zn/Cu: Chemical reaction: \_\_\_\_\_

Direction of e- flow: \_\_\_\_\_

Pb/Cu: Chemical reaction: \_\_\_\_\_

Direction of e- flow: \_\_\_\_\_

Zn/Pb: Chemical reaction: \_\_\_\_\_

Direction of e- flow: \_\_\_\_\_

## B. Concentration Cells

### Calculations:

1. One drop has a volume of about 0.05 mL, thus the concentration of  $\text{Cu}^{2+}$  in the dilute half-cell is about  $(0.05/30) \times 0.1\text{M} = 2 \times 10^{-4}\text{M}$ . Use the Nernst equation to calculate the value of E at 298 K and compare it with your observed value. Why might they differ?

$E_{\text{observed}} =$  \_\_\_\_\_

$E_{\text{calculated}}$  (show calculations) =

---

---

---

---

---

2. After adding 10 drops (0.5 cm<sup>3</sup>) of  $\text{Na}_2\text{S}$  solution ( the concentration of  $\text{S}^{2-}$  is about  $(0.5/30) (0.1\text{M}) = 2 \times 10^{-3}\text{M}$ ), calculate the concentration of  $\text{Cu}^{2+}$  from the relationship of the solubility product:

$$K_{\text{sp}} = [\text{Cu}^{2+}] [\text{S}^{2-}] = 8 \times 10^{-37}.$$

Then calculate the value of E and compare it with the observed value obtained after the addition of  $\text{Na}_2\text{S}$  solution to the dilute  $\text{Cu}^{2+}$  half-cell.

$E_{\text{observed}} =$  \_\_\_\_\_

$E_{\text{calculated}}$  (show calculations) =

---

---

---

## Experiment (11)

### Chemical Thermodynamics


#### Objectives:

- To study the factors that affect the spontaneity of a chemical reaction.
- To determine some thermodynamic functions for a chemical reaction.

#### Introduction:

Several criteria could be used to judge the spontaneity of a chemical change. All these take into account both the enthalpy change of a reaction,  $\Delta H$ , and the change in the entropy or randomness,  $\Delta S$ , of the reaction. One simple measurement is the equilibrium constant,  $K$ .

Consider the reaction:

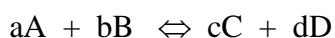

$$K = \frac{[Ce]^c [De]^d}{[Ae]^a [Be]^b}$$

Where  $[Ae]$ ,  $[Be]$ ,  $[Ce]$ ,  $[De]$  are equilibrium concentrations.

At any given concentration, the concentration expression,  $Q$ , is:

$$Q = \frac{[C]^c [D]^d}{[A]^a [B]^b}$$

$K = Q$  only at equilibrium. By comparing  $K$  and  $Q$  values one can predict the direction of the spontaneous reaction. Usually, measurements of  $K$  are not easy to carry out.

A second quantitative criterion of the tendency for reaction to occur is the potential or electromotive force, EMF, of a oxidation-reduction reaction. A positive potential indicates that the reaction will proceed spontaneously in the direction as written and a negative potential indicates that the reaction will proceed spontaneously in the opposite direction.

A third criterion of the tendency of a reaction to occur is the free energy change of the reaction,  $\Delta G^*$ .

Since potentials, equilibrium constants and free energies measure the tendency of a reaction to occur, they must be interrelated, and they are, through the expression:

$$\Delta G^\circ = -nFE^\circ = -RT \ln K$$

Where  $\Delta G^\circ$  is the standard free energy change of the reaction and  $E^\circ$  is the standard cell potential.

If we know  $E^\circ$  for example, both  $K$  and  $\Delta G^\circ$  can be calculated.

Cell potential measurements can often be made quickly and accurately and therefore are favored for finding the tendency of a reaction to occur, especially in aqueous solutions.


In this experiment you will measure EMF and  $\Delta H$  for a chemical reaction. From these you will calculate  $\Delta G$  and  $\Delta S$  values.

---


\* If  $\Delta G$  is negative, the reaction will proceed spontaneously from left to right as written. If  $\Delta G$  is positive, the reaction will proceed spontaneously in the opposite direction.

## Experimental Procedure:

1. Prepare 25 mL of a solution containing 0.10M Fe(NO<sub>3</sub>)<sub>3</sub> and 0.10M FeSO<sub>4</sub>, by mixing 8.50 mL of 0.30M FeSO<sub>4</sub>, 8.50 mL of 0.30 M Fe(NO<sub>3</sub>)<sub>3</sub> and 8.00 mL of H<sub>2</sub>O (using a graduated cylinder).
2. Introduce the 25 mL of this solution into a 150 mL beaker.
3. In another beaker, add about 0.5g of manganese dioxide, MnO<sub>2</sub>, and introduce 25 mL of a ready made mixture of 0.10 M KMnO<sub>4</sub> and 0.1M HNO<sub>3</sub>.
4. Prepare the salt bridge as you did in the previous experiment.
5. Insert the carbon electrodes in the beakers and measure the voltage of the cell consisting of the two half cells:

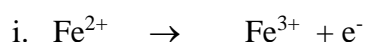

and

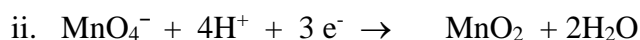

6. When the voltage is stabilized record it.
7. Calculate  $\Delta G$  for the cell reaction using the relationship:
8.  $\Delta G \text{ (kJ)} = -96.5 \text{ nE}$
9. where E is the observed voltage in volts and n is the number of electrons exchange in the overall redox reaction..
10. Measure  $\Delta H$  for the reaction as follows:
  - a) In a 100 mL beaker, place 20mL of 0.20M KMnO<sub>4</sub>, 5.0mL of 6M HNO<sub>3</sub> and 5.0mL of water, (all measured by graduated cylinder).
  - b) Place 20mL of 0.50M FeSO<sub>4</sub> in a flask.
  - c) Carefully measure the temperature of each solution (KMnO<sub>4</sub> and FeSO<sub>4</sub> solutions), to the nearest 0.2°C.
  - d) Add the FeSO<sub>4</sub> solution to the KMnO<sub>4</sub> solution in the beaker with stirring and record the highest observed temperature.

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

**Experiment (11)**  
**Chemical Thermodynamics**  
**Report Sheet**

**Data and Calculations:**

1. Voltage of the cell (V) \_\_\_\_\_
2.  $\Delta G$  for the cell reaction (kJ) \_\_\_\_\_
3. Average temperature of solutions before mixing \_\_\_\_\_
4. Maximum temperature after mixing \_\_\_\_\_
5. Change in temperature ( $\Delta T$ ) \_\_\_\_\_
6. Total volume of final solution \_\_\_\_\_
7. Mass of solution (g) (assume density = 1.0 g/cm<sup>3</sup>) \_\_\_\_\_
8. Specific heat of solution 4.18J/ g°C
9. Heat absorbed by solution (kJ) \_\_\_\_\_
10. Heat absorbed by the beaker (kJ) \_\_\_\_\_  
N.B. Heat capacity of the beaker (100 mL beaker) 30 J/°C  
(150 mL beaker) 40 J/°C
11. Heat evolved by the reaction (kJ) \_\_\_\_\_
12. Calculate  $\Delta H$  per mole of  $\text{MnO}_4^-$  reacted (calculate the initial moles of  $\text{MnO}_4^-$  and  $\text{Fe}^{2+}$  and note carefully which reagent limits the extent of the reaction).  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_
13. From the  $\Delta G$  and  $\Delta H$  values obtained, calculate  $\Delta S$  for reaction under the conditions you used.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## **Experiment (12)**

### **pH, Hydrolysis**

#### **PreLab**

1. Which ion hydrolyzes more,  $\text{Fe}^{3+}$ , or  $\text{Fe}^{2+}$ ? \_\_\_\_\_ . Explain.

---

---

---

2. Write an equation representing the hydrolysis of each ion:

a.  $\text{F}^-$ : \_\_\_\_\_  
b.  $\text{Fe}^{3+}$ : \_\_\_\_\_  
c.  $\text{HSO}_4^-$ : \_\_\_\_\_

3. Predict whether an aqueous solution of each of the following salts is acidic, basic, or neutral:

Salt	Acidic, Basic, or Neutral	Equation justifying your prediction
$\text{Na}_2\text{CO}_3$		
$\text{KNO}_3$		
$\text{AlCl}_3$		
$\text{NH}_4\text{NO}_3$		

## Experiment (12)

### pH, Hydrolysis

#### Objectives:

- To measure the pH of laboratory and home acids and bases.
- To measure the extent of hydrolysis of various ions.

#### Introduction:

##### pH:

The acidity or basicity of most aqueous solutions is frequently due to small concentration of  $H^+$  or  $OH^-$ . The concentration of  $H^+$  or  $OH^-$  is expressed in terms of the pH, such that:

$$pH = -\log [H^+] \quad (1)$$

Water dissociates very slightly producing equal concentrations of  $H^+$  and  $OH^-$  :



At 25°C, the  $[H^+] = [OH^-] = 1.0 \times 10^{-7}$  mol/L producing a pH = 7. Addition of an acid, which either releases  $H^+$  ions or consumes  $OH^-$  ions, increases  $[H^+]$  and thus pH becomes less than 7. Addition of a base, which either releases  $OH^-$  ions or consumes  $H^+$  ions, decreases the  $[H^+]$  in solution, producing a pH greater than 7.

In this experiment, the pH's of some acids and bases will be measured using acid-base indicators or a pH-meter.

##### Hydrolysis:

Salts ionize in aqueous solution forming cations and anions, which are surrounded by polar water molecules (hydrated) to varying degrees depending on the size and the charge of the ion. A strongly hydrated cation attracts the *negative* end of the polar  $H_2O$  molecule and weakens the O-H bond. This liberates free  $H^+$ , producing an acidic solution. The pH is then less than 7. A strongly hydrated anion, on the other hand, attracts the *positive* end of the  $H_2O$  molecule, again weakening the O-H bond and setting  $OH^-$  free, producing a basic solution. The pH is then greater than 7. This splitting of water molecules by ions is called *hydrolysis*.



## Experimental Procedure:

### A - pH measurement with indicators:

1. Clean nine test tubes with soap and water and rinse with distilled water.
2. Place 2 mL of each solution listed in the Report Sheet into separate test tubes and add 5 drops of universal indicator to each. Estimate the pH of each solution.

### pH measurement with pH-meter:


1. Calibrate the pH meter. (This step is done by your laboratory technician).
2. In a 50-mL beaker introduce about 30 mL of each test solution listed in the Report Sheet.
3. Stir the sample, insert the electrodes, and record the pH. Allow sufficient time for the reading to reach a stable value.
4. Remove the electrode and the sample. Place a "waste" beaker under the electrodes, rinse them with distilled water, touch dry with a soft paper towel, and repeat steps 2 and 3 with the other samples described in Part A. Save the test samples for other students to use.

### B - Hydrolysis of Salts:

1. Transfer 30 mL of each solution listed in the Report Sheet into a 50-mL beaker.
2. Measure the pH of each solution with the pH-meter. Identify the ion that hydrolyzes and write an equation for its reaction with water.

### C - Atmospheric Hydrolysis:

Ammonium salts readily hydrolyze with water vapor releasing ammonia gas. The extent of hydrolysis may be detected by smelling the solid crystals. Cautiously smell solid  $\text{NH}_4\text{Cl}$ ,  $\text{NH}_4\text{C}_2\text{H}_3\text{O}_2$ , and  $(\text{NH}_4)_2\text{CO}_3$ . The salt producing the strongest ammonia odor causes the equilibrium

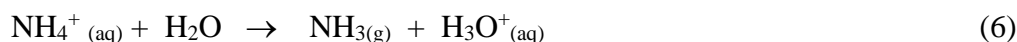

to shift farthest to the right. Therefore, the anion having the strongest affinity for hydronium ion (the strongest base) causes equilibrium (6) to shift farthest to the right, releasing more ammonia gas. Which anion,  $\text{Cl}^-$ ,  $\text{C}_2\text{H}_3\text{O}_2^-$ , or  $\text{CO}_3^{2-}$ , is the strongest base?

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
 ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
 Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (12)

### pH, Hydrolysis

### Report Sheet

#### A. pH

Solution	pH	
	Indicator	pH meter
0.10 M HCl		
0.10 M HC <sub>2</sub> H <sub>3</sub> O <sub>2</sub>		
0.10 M NH <sub>3</sub>		
0.10 M NaOH		
Distilled Water		
Tap Water		
Orange Juice		
Milk		
7-UP		

#### B. Hydrolysis of Salts

Solution	pH	Ion Hydrolyzed	Balanced Equation
	PH-meter		
0.1 M KCl			
0.1 M NaC <sub>2</sub> H <sub>3</sub> O <sub>2</sub>			
0.1 M NH <sub>4</sub> Cl			

**C. Atmospheric Hydrolysis**

1. List the ammonium salts in order of decreasing strength of ammonia smell.

---

2. List the anions in order of decreasing base strength.

---

3. Write an equation for each anion hydrolysis.

a. \_\_\_\_\_

b. \_\_\_\_\_

c. \_\_\_\_\_

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

**Experiment (13)**  
**Acid-Base Titration**

**Prelab**

1-Define the following terms:-

a) Neutralization:

---

---

---

b) Indicator:

---

---

---

c) End Point:

---

---

---

d) Equivalence point:

---

---

---

2-In a titration of NaOH against HCl, 26.12 mL of 0.1624 M NaOH was used to neutralize a certain volume 0.1438 M HCl. What volume of HCl was neutralized?

---

---

---

---

---

## Experiment (13)

### Acid-Base Titration


#### Objectives:

- To determine the standard concentration of sodium hydroxide solution using this standard solution of potassium hydrogen phthalate.
- To determine the concentration of an unknown hydrochloric acid (HCl) solution, using a standard solution of sodium hydroxide.
- To make accurate volume measurements using a pipette, burette and volumetric flask.

#### Notes for good laboratory practice:

1. Always label beakers that contain your solutions.
2. Always rinse out the burette with soap & tap water, distilled water then three times with a small amount of titrant (in this experiment, sodium hydroxide).
3. Do **not** try to get the meniscus to a particular position, e.g. 0.00 mL. A piece of white paper with a black bar placed behind the glass may help you to read the volume.
4. After filling the burette, always remove the funnel. The funnel should **NOT** stay permanently at the top of the burette since any drops of liquid from the funnel will change your readings.
5. Check to see that no bubbles are in the tip of the burette.
6. When titrating, the tip of the burette must be inside the flask as shown below.
7. After draining titrant out of the burette, wait about 10 s and then read the burette.
8. Keep your eye level with the meniscus for all readings.
9. All titrations must be carried out a minimum of three times. The first titration gives an approximate value.

Your hands should hold the stopcock and flask as shown below

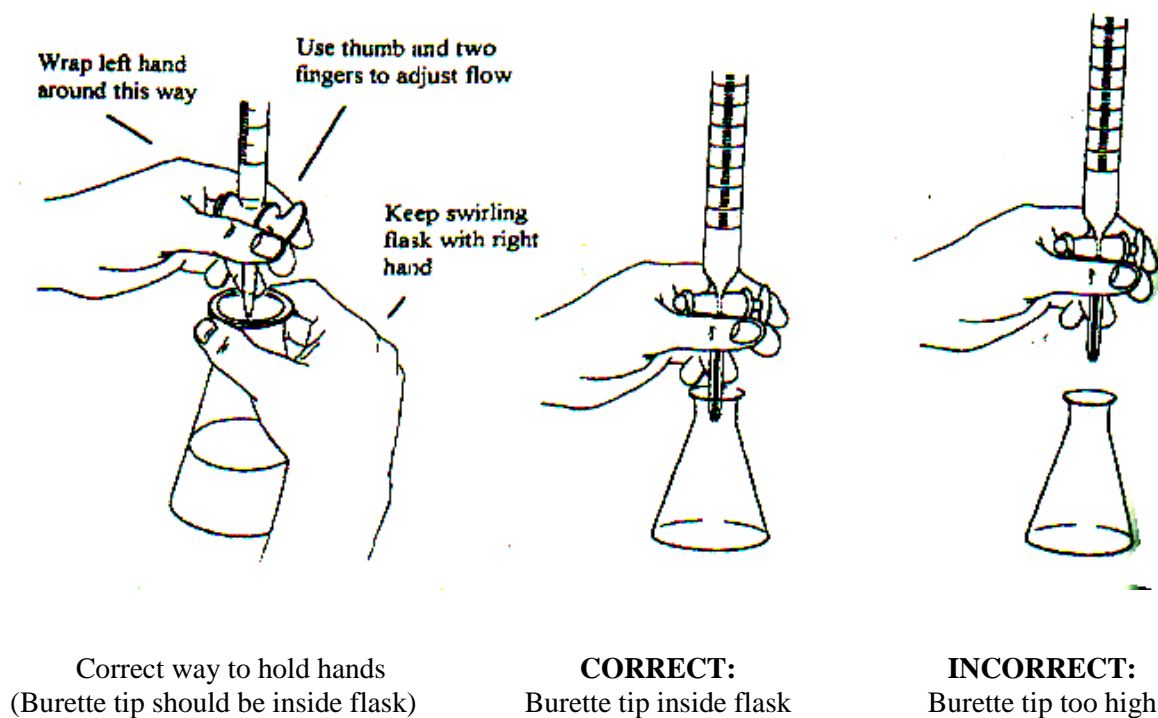

### Experimental Procedure:

#### A. Making the KHP solution:

1. Weigh out, to 2 decimal places, about 0.50 g KHP. Use the spatula to transfer the crystals from the bottle to the 100 mL conical flask. Record the mass in the "Report sheet".
2. Add amount of distilled water to the KHP. Swirl carefully until all the solid has dissolved.

#### B. (1) Standardization of sodium hydroxide solution by titration:

1. Rinse out the burette with a small amount of sodium hydroxide solution. Throw the rinse solution away in the sink.
2. Fill the burette with the NaOH solution, using a small, funnel. Run out enough liquid so that the meniscus is below the zero mark. Make sure that no air bubbles are trapped in the tip of the burette. Remove the funnel from the burette.
3. Record the volume to **two decimal places** as the "Initial burette reading"
4. Add a few drops (about 2 drops) of the phenolphthalein indicator to the KHP solution.
5. Rough titration: Add the NaOH solution from the burette rapidly to the conical flask containing the KHP and indicator. Swirl the flask gently as you add the NaOH so that the solutions mix. When the pink color takes a longer time to disappear, add the NaOH more slowly. **Stop as soon as you see the first permanent light pink color**

6. Record the volume in the burette as the "Final burette reading". Calculate the "Volume added" by subtracting the initial burette reading from the final burette reading. This is the **approximate** volume of NaOH you need for the titration.
7. Rinse out the flask after each titration with distilled water. **Never rinse the flask with KHP solution.** First titration: Refill(if necessary) the burette with NaOH solution and make a note of the new "Initial burette reading" .Prepare your sample as in Steps 1&2 in part A. Take the volume from your rough titration and subtract 2 mL. You can add this volume of sodium hydroxide quickly to your flask. Then add the NaOH slowly and carefully (drop by drop) until the first permanent light pink color is seen. Record the "final burette reading" and calculate the "Volume added".
8. Repeat the titration for a second trial.

### **B. (2) Titration of an unknown (HCl) solution**

**Note:** Ask your instructor for your unknown (HCl) solution

**Repeat** the steps in part B-1 using **unknown ( HCl) solution** instead of the **KHP solution**.

**Repeat** the titration two more times

Finally, **clean all used glassware (including the burette) first with tap water, then with distilled water.**

Student Name: \_\_\_\_\_ Section No.: \_\_\_\_\_  
ID Number: \_\_\_\_\_ Date : \_\_\_\_\_  
Partner Name: \_\_\_\_\_ Group No.: \_\_\_\_\_

## Experiment (13)

### Acid-Base Titration

#### Report Sheet

**Data:**

**Part (B.1): Standardization of a sodium hydroxide by titration:**

	<b>Rough titration</b>	<b>First titration</b>	<b>Second titration</b>
<b>Mass of KHP</b>			
<b>Initial burette reading:</b>			
<b>Final burette reading:</b>			
<b>Volume added:</b>			

Average volume of NaOH solution added = \_\_\_\_\_ mL

**Calculations:**

1. Write a balanced equation for this reaction. **NOTE: The  $\text{KHC}_8\text{H}_4\text{O}_4$  has one acidic proton.**

---

2. Calculate number of moles of KHP.

---

---

---



3. Calculate the number of moles of NaOH.

---

4. Calculate the concentration of the NaOH solution in mole L<sup>-1</sup> in the average volume of the solution used in the titration.

---

---

---

---

**Conclusion:**

The concentration of NaOH solution is \_\_\_\_\_ mole L<sup>-1</sup>. (*To be used in part B.2*)

**Part (B.2) :Titration of an unknown ( HCl) solution:**

	<b>Rough titration</b>	<b>First titration</b>	<b>Second titration</b>
<b>Final burette reading:</b>			
<b>Initial burette reading:</b>			
<b>Volume added:</b>			

Average volume of sodium hydroxide added = \_\_\_\_\_ mL

Number of the unknown hydrochloric acid solution: \_\_\_\_\_

**Calculations:**

1. Write a balanced equation for this reaction.

---

2. Calculate the number of moles of sodium hydroxide in the average volume.

---

---

---

---

3. Calculate the number of moles of HCl in 10.00 mL of the hydrochloric acid solution.

---

---

4. Calculate the concentration of the HCl solution.

---

---

---

---

**Conclusion:**

The concentration of Unknown # \_\_\_\_\_ HCl solution is \_\_\_\_\_ mole L<sup>-1</sup>.

**Experiment (13)**  
**Acid-Base Titration**

**PostLab**

1. What is the effect of the following (increase, decrease or no effect) on the determination of the concentration of HCl?

a) Air bubbles in the burette.

---

---

b) More than the required drops of Phenol naphthalene is added to the flask.

---

---

c) The color at the end point is too pink.

---

---

